Africa

RS Components SA

P.O. Box 12182,
Vorna Valley, 1686
20 Indianapolis Street,
Kyalami Business Park,
Kyalami, Midrand
South Africa
www.rs-components.com

Asia

RS Components Ltd.

Suite 1601, Level 16, Tower 1,
Kowloon Commerce Centre,
51 Kwai Cheong Road,
Kwai Chung, Hong Kong
www.rs-components.com

China

RS Components Ltd.

Suite 23 A-C
East Sea Business Centre
Phase 2
No. 618 Yan'an Eastern Road
Shanghai, 200001
China
www.rs-components.com

Shenzhen Everbest Machinery Industry Co., Ltd.

19th Building, 5th Region, Baiwangxin Industry Park,
Songbai Road, Baimang, Xili, Nanshan, Shenzhen,
China P.C. 518108
Tel: +86.755.27353188
www.cem-instruments.com

Lotus NL B.V.

Koningin Julianaplein 10, 1e Verd, 2595AA,
The Hague, Netherlands.
Tel: +31644168999
E-mail: peter@lotusnl.com

Europe

RS Components Ltd.

PO Box 99, Corby,
Northants.
NN17 9RS
United Kingdom
www.rs-components.com

Japan

RS Components Ltd.

West Tower (12th Floor),
Yokohama Business Park,
134 Godocho, Hodogaya,
Yokohama, Kanagawa 240-0005
Japan
www.rs-components.com

U.S.A

Allied Electronics

7151 Jack Newell Blvd. S.
Fort Worth,
Texas 76118
U.S.A.
www.alliedelec.com

South America

RS Componentes Limitada

Av. Pdte. Eduardo Frei M. 6001-71
Centro Empresas El Cortijo
Conchali, Santiago, Chile
www.rs-components.com

Quick Start Guide

RS-8806S

Stock No: 2010210

Non-Contact Forehead IR Thermometer

1. Features

- Precise non-contact measurements
- User selectable °C or °F
- Selectable Body and Surface temp
- Set Alarm value
- Memorization of the last 32 measurements
- Automatic Data Hold & Auto power off
- Display Resolution 0.1°C(0.1°F)
- Backlight LCD display

2. Normal Temperatures According To Measurement Method

Measurement Method	Normal Temp °C	Normal Temp °F
Rectal	36.6 to 38	97.8 to 100.4
Oral	35.5 to 37.5	95.9 to 99.5
Axillary	34.7 to 37.3	94.4 to 99.1
Ear	35.8 to 38	96.4 to 100.4


3. Normal Temperatures According To Age

Age	Temp °C	Temp °F
0-2 years	36.4 to 38.0	97.5 to 100.4
3-10 years	36.1 to 37.8	97.0 to 100.0
11-65 years	35.9 to 37.6	96.6 to 99.7
> 65 years	35.8 to 37.5	96.4 to 99.5


4. Configuration

- 1-IR Sensor
- 2-LCD Display
- 3-Mode Selection
- 4-Down Button
- 5-Up Button
- 6-Mode Button
- 7-Measurement Trigger
- 8-Battery Cover

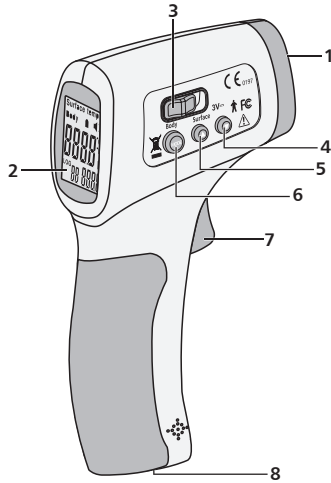5. Indicator

- 1-Surface mode Symbol
- 2-Body mode Symbol
- 3-Digital readout
- 4-Battery Symbol
- 5-The order number
- 6-Save data readout
- 7-Temperature °C(Celsius)/ °F(Fahrenheit) Scale
- 8-Buzzer symbol

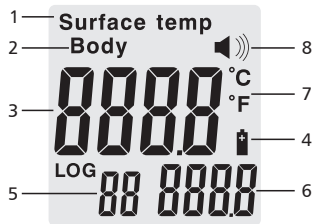6. Technical Specifications

Normal Conditions of Use	
Display Resolution	0.1°C (0.1°F)
Operating Temperature	10 to 40°C (50 to 104°F)
Storage Temperature	0 to 50°C (32 to 122°F)
Humidity Rate	≤ 85%
Power	DC 3V (2 x "AA" batteries)
Size	128 x 74x 36 mm / 5 x 2.9 x 1.4 in (L x W x H)
Weight	Gross 125.4g / Net 104.5g

Measuring Range	
In Body Mode	32.0 to 42.5°C (90 to 108°F)
In Surface Temp Mod	0 to 60°C (32 to 140°F)
Accuracy	±0.3°C (0.54°F)
Emissivity	Fixed at 0.99
Measuring Distance	1 cm – 10 cm (0.39 in – 3.9 in)
Automatic Stop	7 sec.

Non-contact Body Infrared Thermometer Precision

32 to 35.9°C / 93.2 to 96.6°F	±0.3°C / 0.5°F	According to ASTM Standard E1965-1998 (2003)
36 to 39°C / 96.8 to 102.2°F	±0.2°C / 0.4°F	
39 to 42.5°C / 102.2 to 108.5°F	±0.3°C / 0.5°F	

7. Calibration Instructions

For stable and reliable results, The thermometer provides a user-correction function, the operation of the method is as follows:

1st step: take the temperature of a person using a conventional thermometer; you will get 37.5°C (99.5°F) for instance.

2nd step: take the temperature of the same person using the Noncontact Body Infrared Thermometer keeping the 5 to 15 cm (2 in to 5.9 in) distance between the thermometer and the forehead (Take care to remove any obstacle which could alter the measurement (hair, perspiration...)). If you get 37.5°C (99.5°F), the Non-contact Body Infrared Thermometer is properly set and ready for use.

If you get a lower temperature, 36.4°C (97.4°F) for example, your difference is 1.1°C (2.2°F). You should adjust the temperature on the Non-contact Body Infrared Thermometer and add the difference, i.e. 1.1°C (2.2°F).

To do it, press the "MODE" button for 2 seconds, the screen displays F1, press "MODE" button again until you get F3. Press "UP" button in order to add the difference (in our example, 1.1°C – 2.2°F).

3rd step: To check, take the temperature again using the Non-contact Body Infrared Thermometer.


8. Use

8-1. Install battery


8-2. For the first use or when inserting new batteries wait between 10 minutes for the warm-up of the apparatus and when inserting the new batteries.

8-3. If the device is not used for a long time, once you turn it on again, the device will test the room temperature first and will delay turning on for one or two seconds.

8-4. Aim towards the forehead (see the diagram below for the positioning), from a distance of 5cm(2in), press the measuring key, the temperature is displayed immediately. Making sure there is no hair, perspiration, cosmetic or cap covered on the forehead.

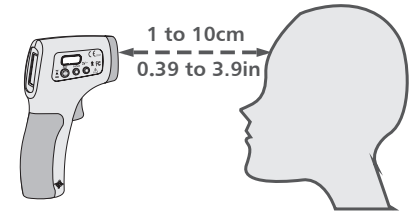8-5. When the room temperature is significant difference, or there is perspiration on the forehead, You can take the temperature behind the ear lobe. Making sure there is no hair, perspiration, cosmetic or cap covered.

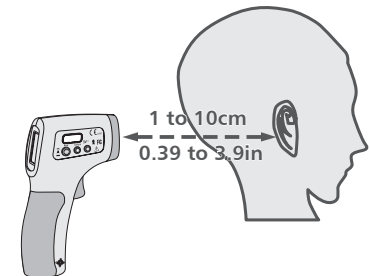9. Measuring Operation

9-1. Choosing The Temperature Unit – F1 Function

Press “MODE” button for 2 seconds, the screen displays: F1. Select “DOWN” for degrees Celsius, “UP” for degrees Fahrenheit.

9-2. Alarm Setup – F2 Menu

Press “MODE” touch for 2 seconds, the screen displays: F1. Press twice “MODE” touch to get F2. Select “UP” to increase the threshold by 0.1°C (0.1°F), “DOWN” to reduce it by 0.1°C (0.1°F).

Note: The alarm threshold default value is 38°C (100.4°F).

9-3. Total Difference – F3 Menu

To adjust the total variation of your the Non-contact Body Infrared Thermometer. Press “MODE” button for 2 seconds, the screen displays: F1. Press “MODE” button two times to get F3. Select “UP” to increase the difference by 0.1°C (0.1°F), “DOWN” to reduce it by 0.1°C (0.1°F). In the cases of seasonal or environmental changes a verification and adjustment should be carried out.

Note: This function is only effective Body

9-4. Buzzer ON/OFF – F4(F3) Menu

Press “MODE” button for 2 seconds, the screen displays: F1. Press three(two) times “MODE” button to get F4(F3).

Select “UP” to open the buzzer (a sound icon “🔊”) is displayed on the LCD screen), press “DOWN” to stop it (the icon disappear).

9-5. Exiting The Setting Mode

Press “MODE” button until the screen turns off.

9-6. The Non-contact Body Infrared Thermometer is specially designed to take the body temperature of a human being. For this, use the Body mode. Measurement range for Body mode: 32 to 42.5°C (86 to 108°F)

You can also use the Non-contact Body Infrared Thermometer to measure the temperature of an area or an object, a food, a liquid or a room temperature. For this, use the Surface mode. Measurement range for Surface mode: 0 to 60°C (32 to 140°F)

Important: The area temperature differs from the internal body temperature. To obtain the internal temperature always use the Body mode.

Please make sure to select the Body mode for an internal temperature reading and the Surface mode for an external area reading (bottle, bath, room...).

9-7. Data Memory

Data memory automatically after temperature measurements, which will display at the right corner of LCD. Press “UP” or “DOWN” button to display the last temperature measurement.

Under the condition of power off, press “UP” and “DOWN” button together for two seconds, to display the last temperature measurement.

At this mode if change the order number to 0 and press “MODE” button, then delete all memory data.

9-8. Changing The Batteries

Display: when the LCD screen displays “🔋”, the battery is used.

Operation: Open the lid and change the batteries, taking great care with the correct positioning. A mistake with this could cause damage to the apparatus and compromise the guarantee of your Non-contact Body Infrared Thermometer. Never use rechargeable batteries. Use only batteries for single usage.

Remove the battery from the instrument if it is not required for extended periods of time in order to avoid damage to the thermometer resulting from a leaking battery.

9-9. Longevity Use

The Non-contact Body Infrared Thermometer was conceived for an intense and professional use, its longevity is guaranteed for 40000 readings.

10. Advice

- The protective glass over the lens is the most important and fragile part of the thermometer, please take great care of it.
- Do not recharge non rechargeable batteries, do not throw in fire.
- Do not expose the thermometer to sunlight or water.

11. Maintenance and Cleaning

- The Infrared Sensor is the most precise part, must be protected carefully.
- Clean the device with a cotton bud lightly moistened with 70% alcohol.
- Do not clean the device with corrosive detergent.
- Keep the device away from water or other liquid.
- Store the device in a dry environment, and keep it away from dust and direct sunlight.

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