



Lohmann & Rauscher

The Debrisoft® product family

Always a good choice.

When I want
to be sure, I

debrisoft it.



When I want effective wound debridement, I

debrisoft it.



When I need rapid results, I

debrisoft it.



Debrisoft® Lolly

- for debridement of deep wounds and hard-to-reach areas
- secure fixation of the fibre head
- ideal for surgical invasive wounds

Debrisoft® Pad
13 x 20 cm

- removal of debris and exudate from the wound
- the angled fibre tips actively loosen debris and exudate from the wound
- absorption and removal of skin keratoses and skin flakes
- protection of newly formed granulation tissue and epithelial cells


Debrisoft® Pad
10 x 10 cm



Why the original? For these reasons:

Debrisoft is **highly effective**, absorbs and removes debris and exudate.


Absorbs, binds and removes bacteria; removes biofilm.

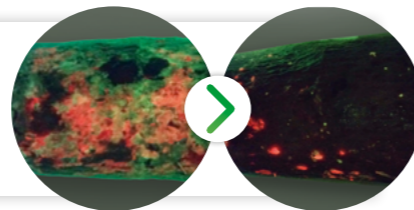

Source 1

Debrisoft is simple to use via its **easy application**.


High user satisfaction; rapidly visible results within 2–4 minutes.

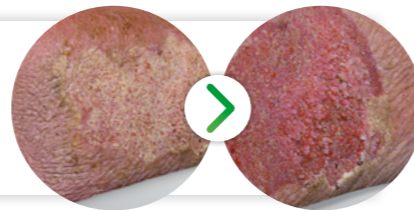

Source 2

Debrisoft means virtually **painless and safe** debridement.

High patient acceptance rate; gentle on intact tissue.

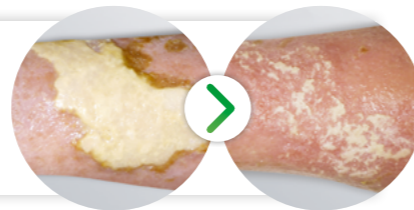

Source 3

First moisten Debrisoft with any common wound cleansing solution (or water).

Debrisoft Lolly's handle is not only ergonomically shaped but also **break-proof**.

Tips & Tricks:

Practical option: To aid moistening Debrisoft, place the irrigation solution onto the transparent packaging layer.

The larger pad is **ideal for burns, 'scaly' legs and large wounds.**



1. Moisten

Debrisoft Pad only cleans effectively if it is moistened beforehand! Debrisoft can be combined with any common wound irrigation solution (or water). Whichever is selected, follow local guidelines.

2. Clean

Applying gentle pressure, clean the wound using the soft fibre side. If necessary, clean the skin surrounding the wound with another moistened Debrisoft Pad.

3. And you are done

Dispose of the used Debrisoft Pad. Always observe the local disposal guidelines.



The Debrisoft® product family

Three options – one solution.



Want to use the original?
Get helpful tips here...



Debrisoft® Pad

sterile, individually sealed

Size (cm)	REF	Shipping units
10 x 10	34 321	5/50
13 x 20	34 323	5/50

Debrisoft® Lolly

sterile, individually sealed

REF	Shipping units
33 224	5/50

Please note the instructions for use.

Experience Debrisoft® for yourself.

Contact Lohmann & Rauscher for your **Debrisoft®** sample pack.

AUS: 1300 572 869 | NZ: 0800 572 869 | Email: info@au.LRmed.com

Lohmann & Rauscher Pty Ltd · 321 Ferntree Gully Road · Mount Waverley, Victoria 3149 · Australia
Lohmann & Rauscher Limited · 21 Queen Street · Auckland 1010 · New Zealand

- 1) Images: Photos used with patient's kind consent, UK.
- 2) Images: Gabriele Stern, ATZ-Stern, Bocholt, Germany.
- 3) Images: Barbara Pritchard, Wrexham Maelor Hospital, Wales, UK.