

REHABILITATION & CLINICAL COUCH RANGE

## **Prioma™** Couch




# Quality by design

#### Safety, comfort and reliability

The Prioma couch range has been designed to provide moving and handling solutions for patients with different mobility levels to assist clinical, medical and therapy staff for rehabilitation and clinical use.

All Prioma couches have been designed to provide clearance for a mobile hoist which ensures safe patient transfer procedures to be carried out.

The range is manufactured to a high standard to ensure safety for both the clinical staff and patients, as well as providing reliability, minimal maintenance and product longevity.

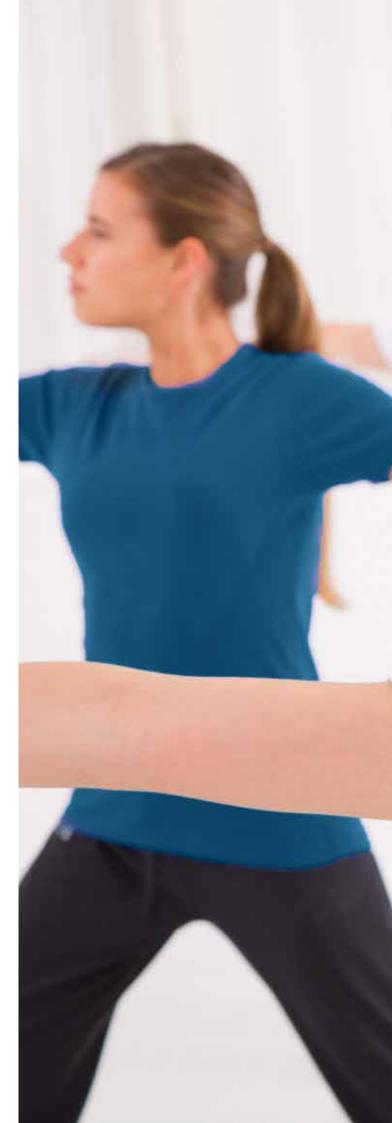



# Prioma electric couch range

#### A fully electric couch offering ergonomic and manual handling solutions

A versatile couch that is suitable for a wide range of examinations, manual therapy techniques and various medical procedures. The versatility of the 2-section and 3-section models allows greater flexibility when positioning patients for a wider range of examination and treatment procedures.

The optional breathing hole with removable plug, allows face down treatment and with a range of optional features and accessories, the Prioma Electric Couch can cater for most Clinical and Rehabilitation needs.



#### Standard features

- Safe Working Load (SWL) 200 Kg
- Electrically operated variable height adjustment
- Electrically operated backrest positioning (Prioma EC2 & Prioma EC3 only)
- Electrically operated footrest positioning (Prioma EC3 only)
- Tente 125mm central locking castors
- ABS crank arm covers
- M1 flame retardant vinyl upholstery
- Compatible with mobile hoist

#### **Optional features**

- Emergency battery back up unit
- Breath hole, located on the backrest (Prioma EC2 and Prioma EC3 only)

#### **Optional accessories**

- Paper roll holder
- Long side support rails pair (Prioma EC1 only)
- Short side support rails pair (Prioma EC2 and Prioma EC3 only)



**Prioma EC1** 




**Prioma EC2** 



#### Prioma EC1


1-section, adjustable height

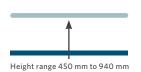




#### **Prioma EC3**


3-section, adjustable height, backrest and foot positioning

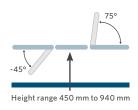


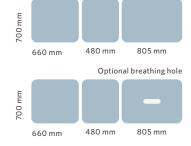

#### **Prioma EC2**

2-section, adjustable height and backrest positioning

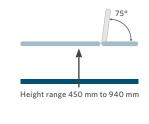

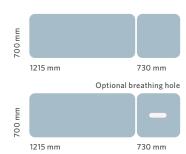





Safe working load (SWL)	200 Kg
Maximum patient weight (MPW)	165 Kg
Accessories weight	
Short side support rails - pair	8 Kg
Long side support rails (Prioma EC1 only) - pair	12 Kg
Paper towel rail	1.5 Kg

# Prioma semi-electric couch range

A semi-electric couch offering ergonomic and manual handling solutions

Available in 1, 2 or 3 sections with electrically operated height adjustment providing ease of use at the touch of a switch.

The 2 and 3-section models come with gas-assisted positioning of the backrest and footrest to cater to a broader range of clinical and rehabilitation needs.



#### Standard features

- Safe Working Load (SWL) 200 Kg
- Electrically operated variable height adjustment
- Manually operated gas assisted backrest (Prioma HC2 and Prioma HC3 only)
- Manually operated gas assisted foot rest positioning (Prioma HC3 only)
- 75mm retractable wheels
- Retractable wheel base
- Levelling foot for uneven floors
- M1 flame retardant vinyl upholstery
- Compatible with mobile hoist

#### **Optional features**

- ABS undercarriage cover
- Breath hole, located on the backrest (Prioma HC2 and Prioma HC3 only)

#### **Optional accessories**

- Paper roll holder
- Long side support rails pair (optional for Prioma HC1 only)
- Short side support rails pair (optional for Prioma HC2 and



**Prioma HC1** 






**Prioma HC2** 

#### **Prioma HC1**


1-section, adjustable height

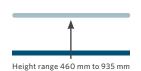




#### **Prioma HC2**


2-section, adjustable height and backrest positioning

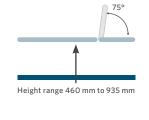


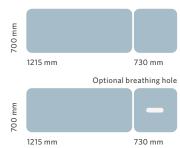

#### **Prioma HC3**

3-section, adjustable height, backrest and foot positioning

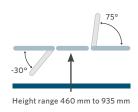

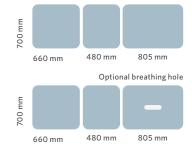





Safe working load (SWL)	200 Kg
Maximum patient weight (MPW)	165 Kg
Accessories weight	
Short side support rails - pair	8 Kg
Long side support rails (Prioma HC1 only) - pair	12 Kg
Paper towel rail	1.5 Kg

# Feature descriptions



### Retractable wheel base and levelling foot

The retractable wheel base ensures ease of transport when the couch needs to be moved. The adjustable foot provides stability when being used on an uneven surface.



### Linak hi-low drive with footswitch

Hands-free control, designed for ease of use. Allows therapy staff to select the correct positioning of the couch, whilst attending to the required moving and handling needs of the patient.



Gas assisted positioning

Prioma semi-electric couch models only.

Gas spring release handles located
on either side of the couch for manual
positioning of the backrest and/or foot end.

### Standard colour options

The vinyl is suitable for a wide range of healthcare applications. The surface is hard wearing, waterproof, easy to clean and is available in 3 standard colours.



**Stone Blue** 



Pebble Black



**Cloud Grey** 

Non-standard colours available and will be priced on application.





This symbol denotes the Safe Working Load of the couch. Maximum Patient Weight (MPW) recommendation takes into account the additional use of accessories on the couch, and subtracts this potential weight from the SWL figure to arrive at a MPW recommendation. If the combined weight of patient and accessories exceeds the noted SWL weight for each couch, the maximum patient weight must be reduced accordingly.



All Prioma couches featured in this brochure have been designed to provide clearance for a mobile hoist to support patient handling and transfer procedures.



All height range measurements are from floor to top surface working height and have a tolerance of +/-10 mm. All angular measurements have a tolerance of  $+/-5^{\circ}$ . The foot section angle depends on the height of the couch as it automatically adjusts to avoid contact with the undercarriage. The angles indicated are the maximum range and based on the couch set at a relative height. All sizes for each individual section do not include hinges and have a tolerance of +/-10 mm.



The upholstered couch sections satisfies flammability requirements to the following standards: BS7176:2007 (medium hazard) using ignition sources 0 (cigarette), 1 (gas) as defined in BS EN. 597-1 & 2 and 5 (crib) BS.5852:2006.

#### **Electrical data**

Power in - 220V AC 50/60Hz, Max current 2.5A

Duty rating - intermittent 10% max 2 min / 18 min

Electric Shock Protection - Class 2 Type B

Liquid ingress protection - IPX4

Internally powered equipment - battery backup only - 2 x 12V series connected batteries, sealed, rechargeable lead/acid gel 1.2 amp hrs



At Arjo, we are committed to improving the everyday lives of people affected by reduced mobility and age-related health challenges. With products and solutions that ensure ergonomic patient handling, personal hygiene, disinfection, diagnostics, and the effective prevention of pressure ulcers and venous thromboembolism, we help professionals across care environments to continually raise the standard of safe and dignified care. Everything we do, we do with people in mind. www.arjo.com

Arjo AB, Hans Michelsensgatan 10, SE-211 20 Malmö, Sweden. Phone: +46 (0) 10 335 45 00

www.arjo.com

