

## Give your mobility a boost.

While it looks and operates mostly like a manual wheelchair, BoostGlide has a power assist motor to make movement almost effortless. When the hand rim is pushed, the clever drive motor intuitively takes over much of the work. It makes uphill climbs a breeze, plus there's a power brake for downhill safety. The battery allows journeys of up to 15km on a full charge.

### Completely customisable.

BoostGlide is a fully-scripted wheelchair, built to the user's precise requirements in our Australian workshop.

Almost endless customisation can be specified to provide improved comfort and support, from footplate or armrest options to a wide selection of colour choices.

BoostGlide wheelchairs can also be adjusted by the user to change seat positioning vertically and horizontally, the footplate height can be changed, and the armrests can be lowered, lifted, or removed.

## Easy climbs. Longer journeys.

With a powerful and reliable drive motor to offer assistance, plus an aluminium frame, it's easy to see why BoostGlide is the leading Australian manufactured power assist wheelchair.



### Customised seating.

For ultimate comfort, many seating options available. Choose customised cushion, backrest and headrest.



### Easy folding mechanism.

By lifting the seat upholstery from the centre, the BoostGlide folds flat enough to fit into a standard car boot space.



### Wheels.

BoostGlide has 24" pneumatic rear wheels with precision bearings for a smooth ride.



### Adjustable armrests.

The armrests have a button lock system to adjust the height to user preference.



### Strong, yet light.

The versatile aluminium folding frame is engineered for strength, yet keeps the overall chair weight low.



### Colour choices.

Choose from an extensive colour palette for the ultimate customisation.

## BoostGlide offers day-long power assistance.

### Features.

- Strong lightweight aluminium frame
- Easy folding for transport or storage
- Hard-wearing upholstery
- Rear wheels 24" x 1 3/8" pneumatic
- Front castors 175mm x 40mm MCP PU (solid)
- Adjustable armrests
- Adjustable footrests

### Safety.

- Crash tested to AS 3696:19 to max user weight 125kg

## Arrange a trial.

Chat to a Glide wheelchair specialist and experience the extraordinary features of the BoostGlide in person. Find your nearest distributor at [glide.com.au/nearest](http://glide.com.au/nearest)

## Specifications.




**Maximum User Weight**  
Up to 125kg



**Overall Length**  
102cm



**Seat Widths**  
37 to 52cm

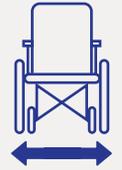

**Overall Width**  
60cm (wheel base)



**Seat Depths**  
33 to 52cm



**Seat Height**  
46 to 53cm



**Chair Weight**  
30.4kg



**Arm Rests**  
Adjustable



**Leg Rests**  
Swing away



**Brakes**  
Electromagnetic  
(via hand rims)  
and manual



**Frame Colour**  
Wide choice



[www.glide.com.au](http://www.glide.com.au)

(08) 9345 3400 | [sales2@glide.com.au](mailto:sales2@glide.com.au)

25 Ledger Road, Balcatta, Western Australia 6021

@GlideProductsAustralia

Glide Products