



IGNITERF



THE INMODE ERA



IGNITERF

The Ultimate non-surgical, minimally invasive, and high-definition skin tightening and contouring platform.

Consumer quest for body contouring treatments continues to be at an all-time high with patients seeking minimal downtime and natural-looking results in their desire to achieve a more youthful appearance. IgniteRF, with cutting-edge technological advancements, provides you with precision control, greater functionality, and more power, allowing you to deliver the enhanced clinical results your patients expect.

IgniteRF is an advanced body contouring and shaping platform, leveraging multiple modalities that synergistically blend to provide enhanced speed and efficacy. With IgniteRF, physicians can deliver customised, combination procedures putting more innovative control in your hands.

Welcome to The InMode Era.

**REFINED TECHNOLOGY.
ENHANCED PRECISION.
EXCEPTIONAL RESULTS.**

Introducing QUANTUMRF

Exceptional skin-tightening and contouring without aspiration.
FDA cleared for coagulation and contraction of soft tissue or where hemostasis is needed.


25

QuantumRF 25 is an advanced technology discharging fractionated radiofrequency (FRF) to maximal depths. The Teflon coated, lightweight handpiece, delivers bipolar pulsed RF to deep body areas without damaging the dermis, while reducing the risk of post inflammatory hyperpigmentation (PIH).



10

QuantumRF 10 is an innovative technology delivering fractionated radiofrequency (FRF) to maximal depths in smaller body zones. The small, easy-to-use handpiece delivers precision bipolar pulsed RF energy to the subdermal plexus, resulting in shorter downtimes and safer procedures.

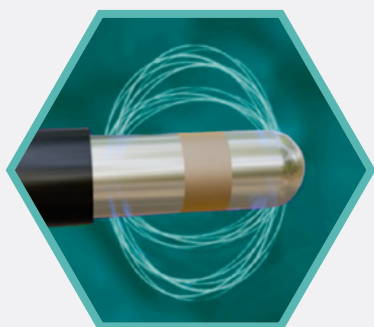

<i>Specifications:</i>	25	10
Cannula Length:	25cm	10cm
Cannula Diameter:	3.3mm	2.2mm
Peak Power:	60W	35W
Maximum Pulse Width:	250ms	250ms

Advantages:

- » No cutting beneath the skin
- » Non-surgical doctor-only procedure
- » In-clinic based procedure
- » Anaesthesia via mesocaine/multi-injector
- » No compression garment needed
- » Low downtime
- » Possibility of therapeutic combination with selected surgical and aesthetic treatments
- » Integrates with a safety mechanism to enhance efficiency and standardise results

LEADERS IN CLINICAL EXCELLENCE



Introducing BODYTITE *turbo*



BodyTite turbo, engineered with enhanced materials, is equipped with a more rigid cannula, ensuring the cannula remains firm when targeting specific body areas. With 50% more power and a 40% longer cannula, BodyTite Turbo delivers even more power, more precision, and faster results.

FACETITE *turbo*

FaceTite Turbo is a non-surgical doctor-only procedure. This applicator is designed with a rigid cannula, providing greater control, precision, and safety in small, delicate tissue zones.



ACCUTITE

AccuTite is the smallest minimally invasive RF device in the Injectable Dermal Radiofrequency family of products. It is an advanced technology for the precise treatment of delicate areas through soft tissue coagulation. Built-in safeguards include real time measurements of skin temperature, impedance monitoring, power cut-off, and audible feedback.

Advantages:



- » No cutting beneath the skin
- » In-clinic based procedure
- » **Single treatment**
- » Long term effect
- » Small amount of anaesthesia or tumescent depending on the treatment area (<100ml)
- » Possibility of therapeutic combination with selected surgical and aesthetic treatments.



IGNITERF - ADVANCED CONTROL & RELIABILITY


- » Full spectrum treatment depths from superficial at 0.5mm to deep treatment of 50mm
- » Preset energy parameters control heat dispersion for optimal treatment results
- » Controlled, even energy distribution delivers uniform large volume heating
- » Temperature feedback at 1,000 times/second ensures immediate real time monitoring
- » 100% energy concentration in the focal treatment zones, protecting peripheral nerves and structures

BODYTITE *turbo*

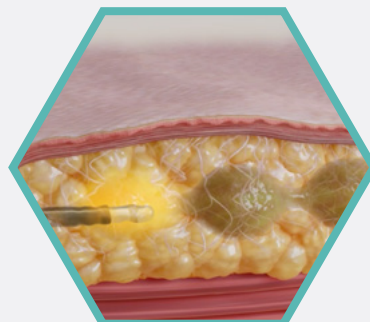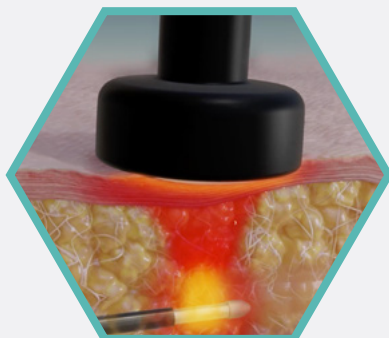

QUANTUMRF

MORPHEUS *burst*

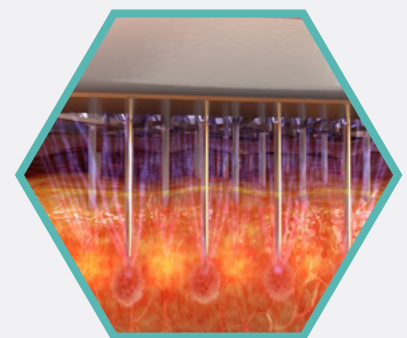

MORPHEUS8 *with*



burst

First and Only FDA cleared Fractional Radiofrequency Microneedling (FRF) technology for Soft Tissue Contraction.


Morpheus8 burst delivers fractional radiofrequency (FRF) energy for clinical applications which require contraction/coagulation of soft tissue or hemostasis. Equipped with proprietary Burst and Scale functionalities, Morpheus8 burst has the unique functionality to automatically deploy bipolar RF to multi-level soft tissue depths in a single pulse, contracting both adipose tissue and dermal tissue. This modular portfolio of tips, including Prime 12 pin, Morpheus8 24 pin, Resurfacing 24 pin and Morpheus8 Deep's large 40 pin, deploys bipolar RF energy in millisecond intervals, allowing you to deliver enhanced clinical results for your patients.

burst deep



Specifications:

	MORPHEUS8 BURST DEEP	MORPHEUS8 BURST	RESURFACING	PRIME
Burst Depth:	multiple depths between 2mm - 7mm	double depth between 2mm - 5mm		
Scale Depth:	multiple depths between 2mm - 7mm	double depth between 2mm - 5mm		
Number of Micropins:	40	24	24	12
Depth Profile:	2mm - 7mm	1mm - 5mm	0.5mm	1mm - 4mm
Heat Profile:	+1.0mm	+1.0mm	+0.25mm	+1.0mm


The Evolution of RF Microneedling

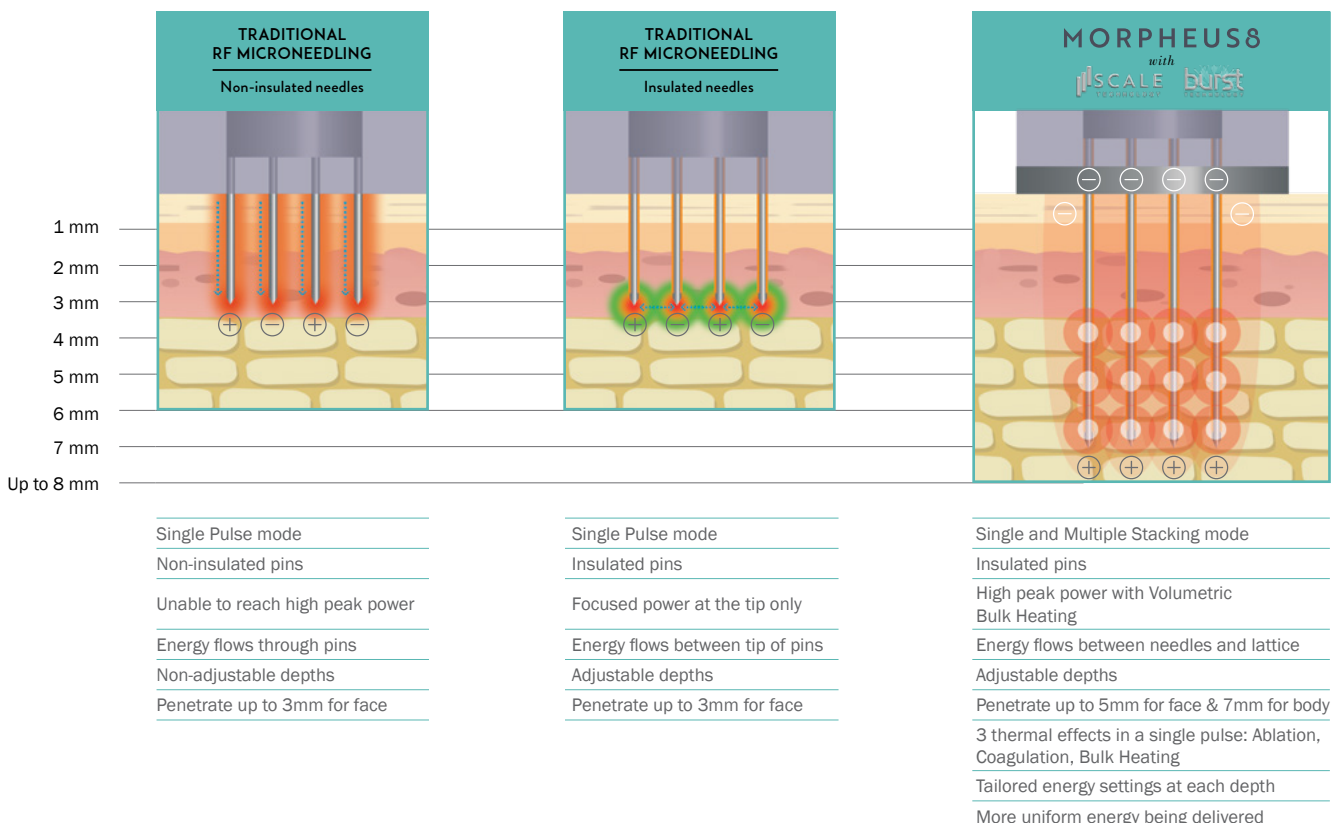


BEFORE & AFTER

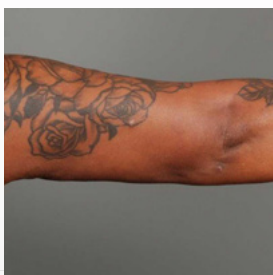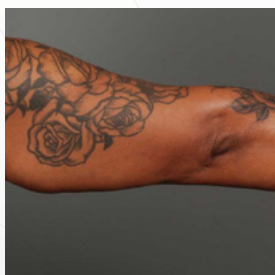



DR. J. RANIERE, JR.

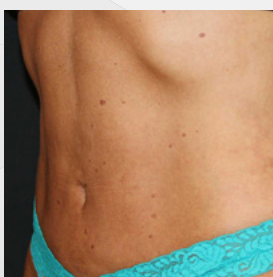

DR. M. LOFFREDO + DR. S. JONES


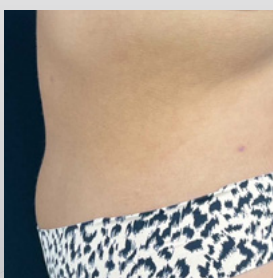


DR. P. HESTER

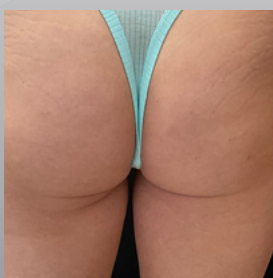
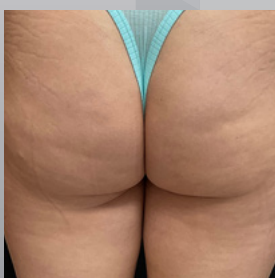


DIAMOND AESTHETICS

BEFORE & AFTER

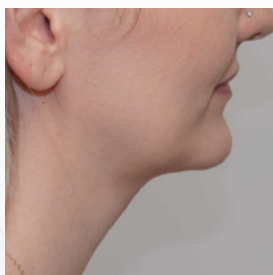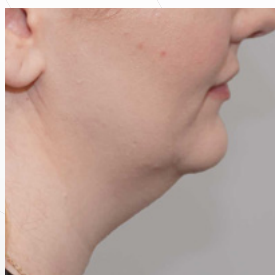

DR. M. TARAJKI

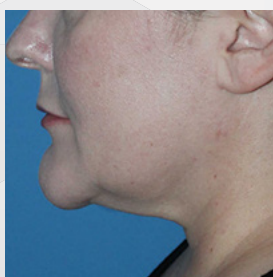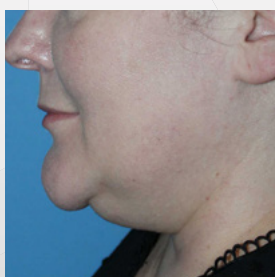

DR. R. DIEPENBROCK






DR. Z. WILLIS



DR. R. MALHOTRA

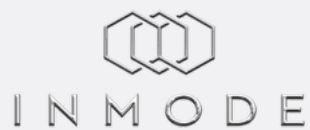REFINED TECHNOLOGY.
ENHANCED PRECISION.
EXCEPTIONAL RESULTS.



 InMode Australia  InMode_Australia  InModeAustralia  盈美特 InMode Australia  458006879

INFO@INMODE.COM.AU | INMODE.COM.AU | +61 (2) 7209 5678