
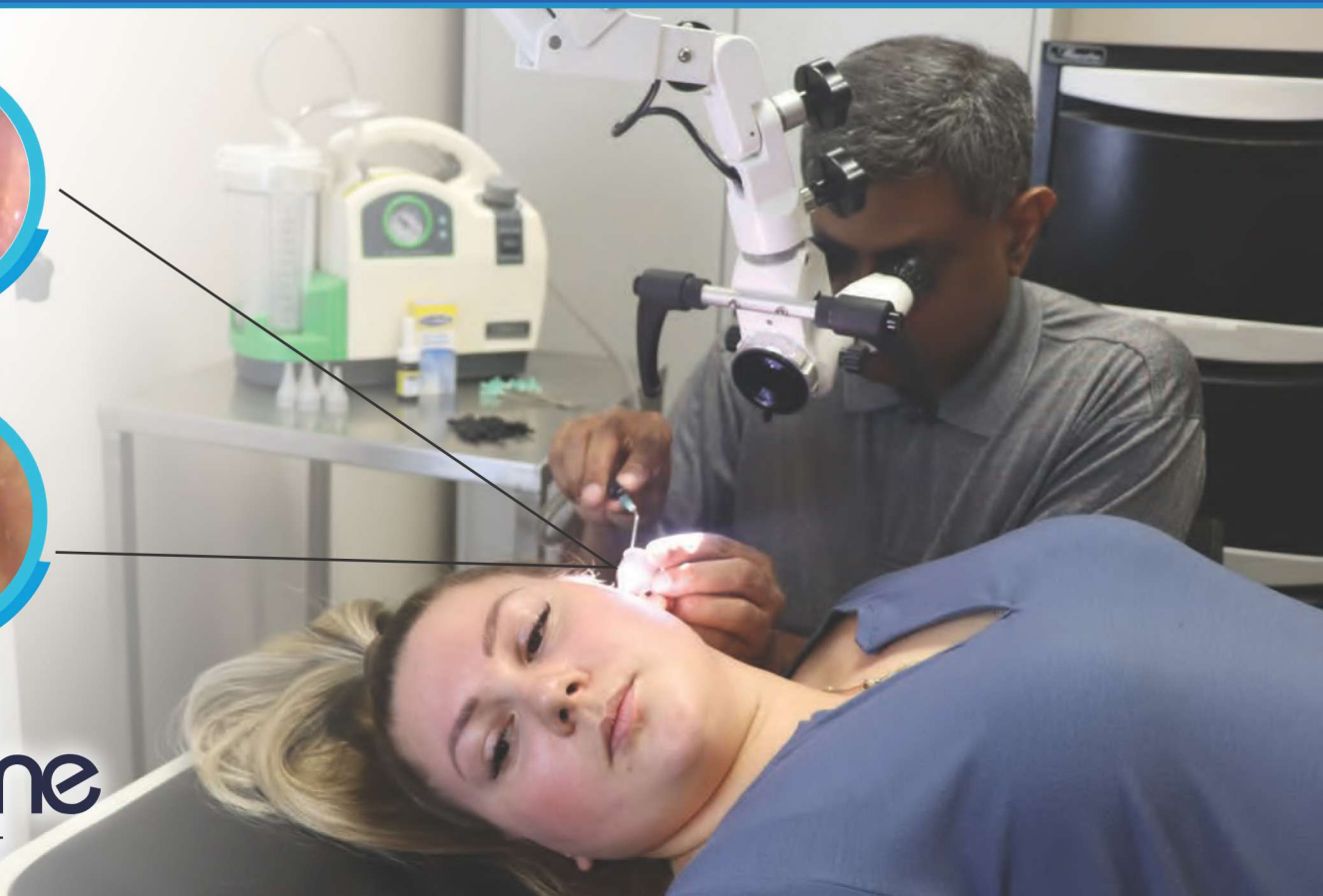
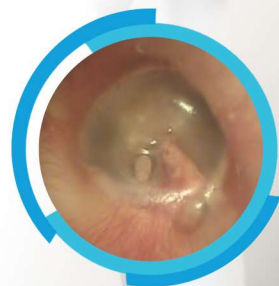


# Ear Toilet Pack

---

## Microscopic Suctioning

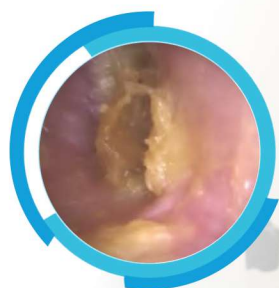

## WHAT IS EAR TOILETING ?

Ear cleaning with micro suction (Ear Toileting) is a progressive treatment which allows the removal of debris and ear wax with a small suction tube attached to a suction controller inserted directly into the ear canal. The thin tube gently vacuums debris out while the whole ear area is under inspection with a microscope by the practitioner.

Micro suction relies on a vacuum movement to remove ear wax or debris and is one of the cleanest, simplest and safest ways to remove ear wax build up and clear blockages in the ear canal.




## CUSTOMER BENEFITS


Micro suction is accepted across the medical profession as the safest and most comfortable method of debris and ear wax removal. With no liquids being used there is no mess during the procedure and won't lead to further damage that other ear cleaning procedures may cause such as ear syringing or ear candling. The fact that the ear canal is being observed throughout the process, makes it exceptionally safe. It is also a comfortable process and is usually undertaken in a few minutes. Most people who have had their ears cleaned of wax using this proven method say it was the gentlest and best ear cleaning experience they have had.

- Fast procedure time, results are instant
- Well tolerate, comfortable for the patient
- Safest method of ear cleaning
- Clean & hygienic

BEFORE

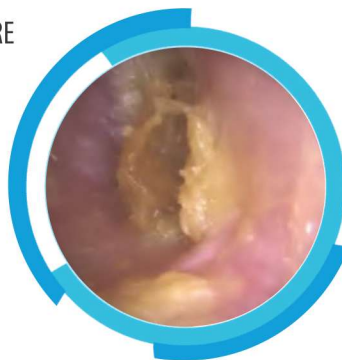


AFTER

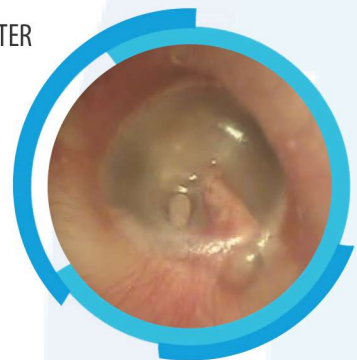

## YSX LED ENT OPERATING MICROSCOPE

Precise structure, modern appearance with a LED illumination system that offers excellent parfocal, large field view, precision and clear three dimensional images. The YSX series of operation microscopes have optional working distances and magnification and are suitable for numerous departments including ENT Specialists, Operating Theatres, Dentists, General Practices, Ophthalmology & Gynaecology

- 1500mm Super long pantograph arm offers a wide operation space
- Four different working distances: 200mm (standard), 250mm, 300mm, 400mm
- Smooth turret 3-step magnification pod 6x, 12x, 20x.
- Field diameter: 7.88mm-81mm
- Magnification: 2.4X-12X
- Orange, Green & Clear Filters
- 10W LED illumination
- 20, 000 hours LED life span
- 50, 000 Lux object illumination

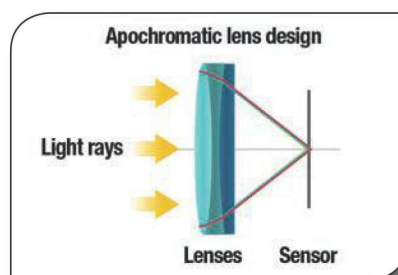

RACGP

*"Other options for mechanical removal include microsuction, with or without curettage, under microscope. Unfortunately, this is usually only available in dedicated ENT outpatient clinics. Microsuctioning is generally well tolerated, safe and efficacious. It also has the advantage of not exposing the ear to moisture and thus has fewer contraindications and is associated with a lower frequency of infections. Discomfort due to noise is the greatest complaint. However, there is no clinical evidence that this affects hearing."*

Visit the RACGP website & search "Ear Wax Management" for more information

# LED YSX FEATURES

Precise spring balanced arm for effortless positioning & comfortable operation

Adjustable Illumination for comfortable visualisation

Uniform series of joints for full 360 degree spherical movement

Premium German made Optics for precise inspection

1500mm Super long pantograph arm offers a wide operation space

Max: 1600mm

3-step magnification pod (6x, 12x, 20x)

Max 1200mm


Steel base for guaranteed stability

Max: 1550mm

Min: 600mm

Smooth lockable castors for fixed secure operation

Base: 600mm

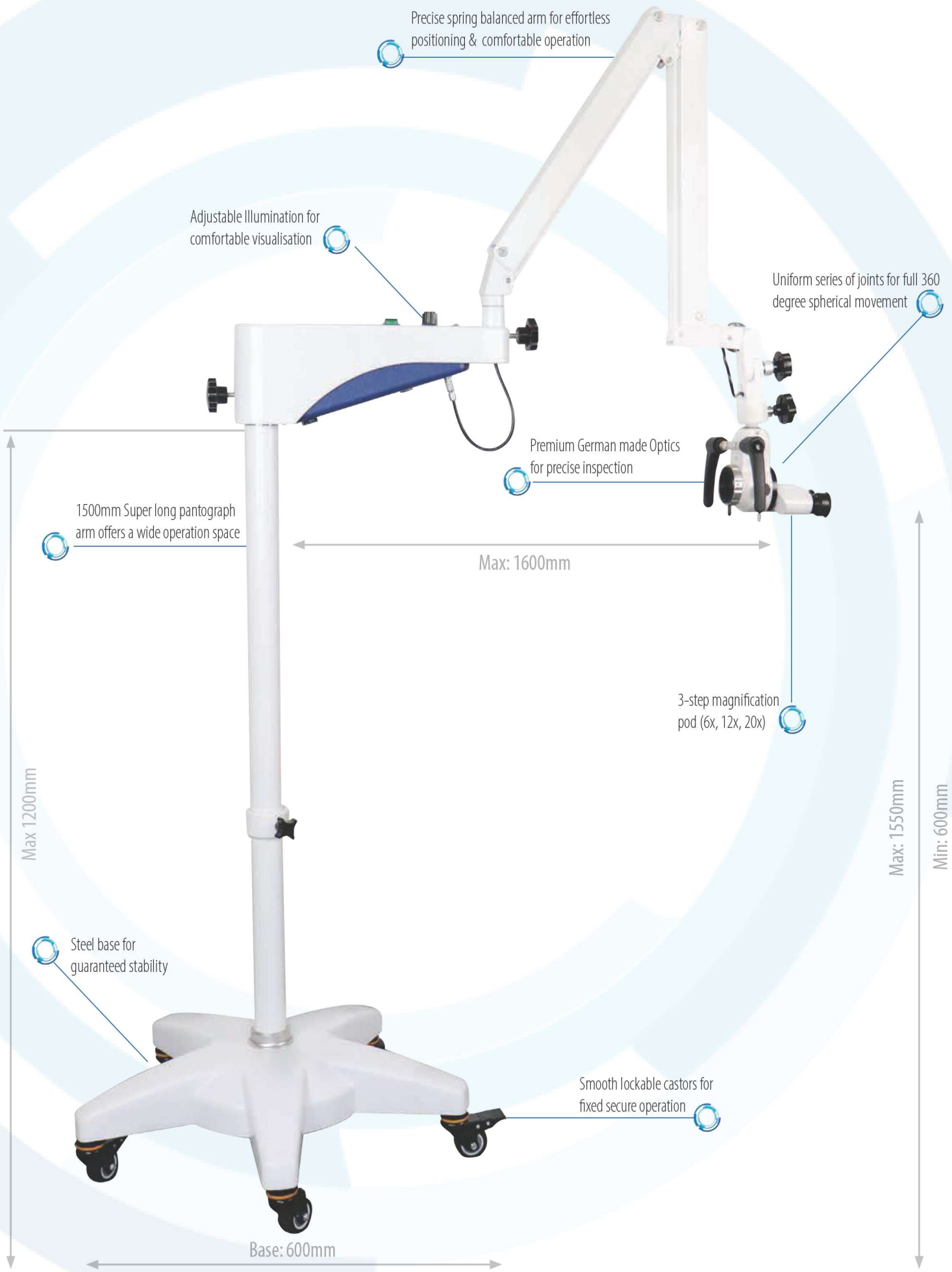

## INCLUSIONS



### TONGYE SUCTION UNIT

Portable suction Machine with Filter, Collection Canister and Aspiration Tubing



### EAR SPECULUM

Precise positioning & easy operation.



### SUCTION TUBES

Ergonomic Micro Suction Tubes. Packs of x 100. Clean. Disposable. 14 Gauge.



### CONTROLLER

Adjust the level of suction through thumb vent mechanism for precise suction control



### WAX REMOVER

Breaks down hard wax to make removal via suction effortless.

**zone**  
MEDICAL

Zone Medical PTY LTD  
[www.zonemedical.com.au](http://www.zonemedical.com.au)  
Phone: 1300 009 663  
Email: [info@zonemedical.com.au](mailto:info@zonemedical.com.au)  
7/22 Mavis Crt Ormeau, QLD 4208