



MedLink™

TouchPoint™
Mobile Technology Cart



Humanscale®
HEALTHCARE



Scan a patient wristband or
Scan a location barcode or
enter a device.

Area	SCOP	Info
1	2	3
4	5	6
7	8	9
10	11	12

DELL

About Humanscale Healthcare

We believe that better design leads to better care.

We use our ergonomic expertise—developed over 30 years of improving efficiency and comfort in the workplace—to create next generation medication administration carts, mobile workstations, wall mounts, seating and accessories for the healthcare environment.

All our solutions are designed with function, simplicity, and longevity in mind. We are experts in designing products that adapt to people's work and their environment.

Our design-led approach frees caregivers to focus on what matters—patient care. Our products actively streamline everyday workflows. They are effortless, easy to use, and make it easy to interact with patients, leaving caregivers free to prioritize their patients, and offer better care.

Better Design, Better Care.



The MedLink Mobile Technology Cart is a point-of-care solution that ensures safe, accurate medication administration. MedLink™ is specifically designed to reduce human error, keep patients safe, and create a closed-loop medication system.

IMPROVING PATIENT SAFETY

Administration errors account for 26% to 32% of total medication errors. By linking medication administration to patient wristbands, MedLink significantly reduces the risk of error during the administration process.

MedLink creates a closed-loop medication system. First, individual drawers are loaded with medication at the pharmacy, ADC, or nursing station. Then each drawer is electronically labeled in the MedLink system and assigned to a specific patient. When that patient's wristband is scanned, the drawer unlocks automatically.

INCREASING ACCOUNTABILITY

Every event is recorded—log-ins, scanning, overrides and drawer opening/closing—introducing a new level of accountability and visibility into the medication workflow.

STREAMLINING WORKFLOWS


MedLink offers significant efficiency gains for any healthcare organization. It streamlines the medication transportation workflow, and reduces the cost of injury claims and medication errors. MedLink's intelligent automation optimizes medication administration, so caregivers can focus on what truly matters—patient care.

Closed-Loop Medication System

MedLink creates a closed loop medication system using lockable drawers and patient wristband scanning. This ensures each patient receives the correct medication, while also allowing caregivers to prioritize patient care over complex administration.

HOW IT WORKS

1. Fill the MedLink cart at a secure, designated filling area.
2. Scan patient wristband to access medication.
3. Patient medication drawer opens automatically, ready for administration.

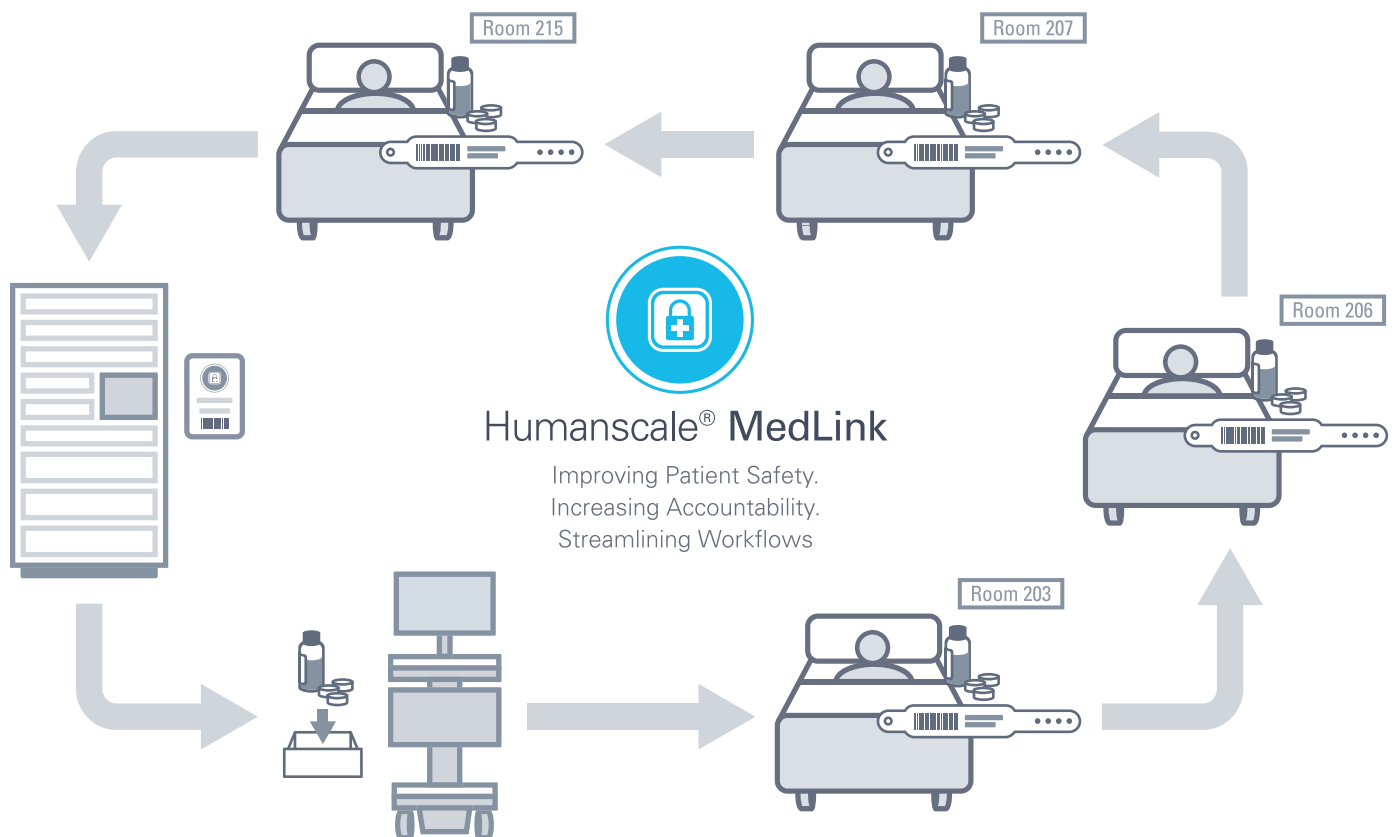

MedLink Hardware

MODULAR. CONFIGURABLE. SIMPLE.

MedLink's modular system can be upgraded and adapted to suit your workflow. MedLink Housings include two removable cassettes that can be transported and exchanged while drawers remain locked and secured.



Standard

M1 is the standard configuration that can support 16 small drawers, 8 medium drawers, 4 large drawers — offering numerous drawer combinations.



High-Capacity

M2, the high-capacity configuration, can support 32 small drawers, 16 medium drawers, 8 large drawers, or a combination of sizes customized to your organizational needs.



Drawers

CUSTOMIZABLE

MedLink's modular drawer design allows customizable drawer layout and sizes to best meet workflow demands.*



Small:
2.9" W x 6.4" D x 2" T




Medium:
6.4" W x 6.4" D x 2" T



Large:
6.4" W x 6.4" D x 4.6" T



Extra Large:
6.4" W x 6.4" D x 7.1" T

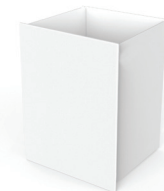

Double Extra Large:
6.4" W x 6.4" D x 9.6" T

DRAWER SIZES

FEATURES

12 different drawer kit configurations available

Notched divider design enables caregivers to lay syringes flat inside drawers while still separating drawer contents

Tool-less removal of shelves, allowing rapid drawer reconfiguration based upon changing workflow requirements


Transparent dividers give caregivers the ability to view medications in rear drawer sections and customize the internal layout of the drawers

RFID embedded technology assists in locating each patient drawer—to the specific cart, cassette and caregiver




*Drawers, shelves and dividers are sold separately, both in kits and individually.

MedLink Lite Software

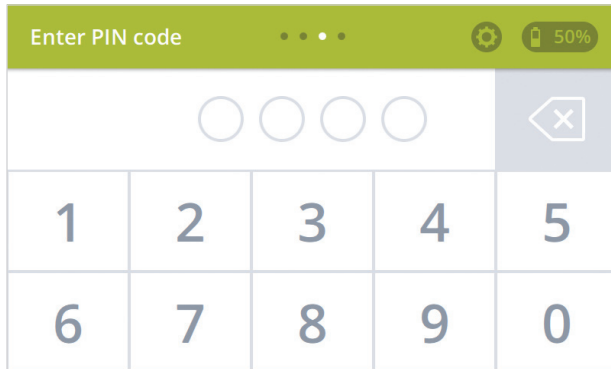


USER INTERFACE AND SECURITY

- Integrated into the User Interface on the T7 dashboard for streamlined control of drawers
- 4-digit PIN code restricts access to drawer operation

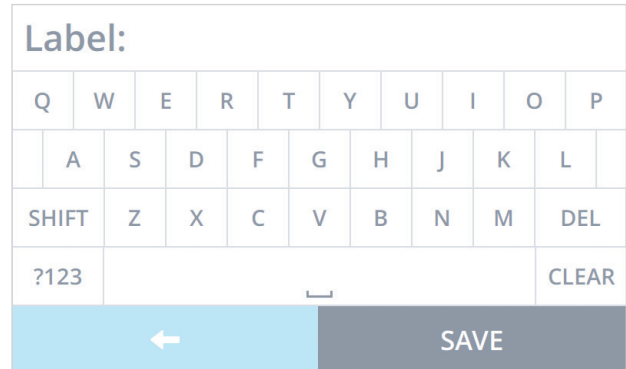

ELECTRONIC DRAWER LABELING

- Patient drawers are labeled through the GUI interface
- Drawers can be labeled by patient name, room number or tagged for supplies

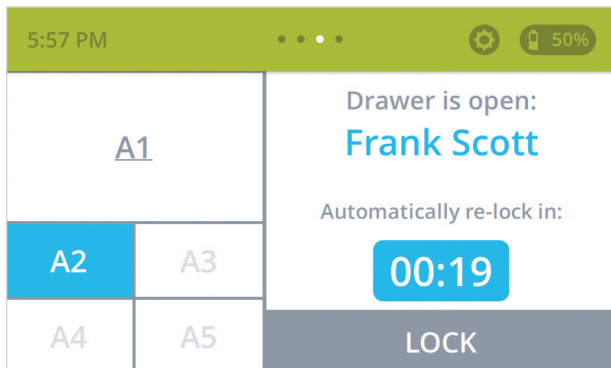


AUTO-CLOSE

- If patient drawers are left open, the system will automatically close and lock the drawer, preventing unauthorized access to medication if the caregiver is called away

MedLink Pro Software

Name	Gender	Date of Birth	Room
Harrison, Holly	F	10/10/1934	RM 108
Harrison, John	M	10/06/1943	RM 104
Jones, Kim	F	10/09/1919	RM 107
Landers, Alex	M	10/15/1978	RM 102
Lang, Amanda	F	10/08/1987	RM 106
Richards, Ricardo	M	10/05/1990	RM 103
Scott, Frank	M	09/15/1956	RM 101
Smith, Sally	F	10/07/1964	RM 105


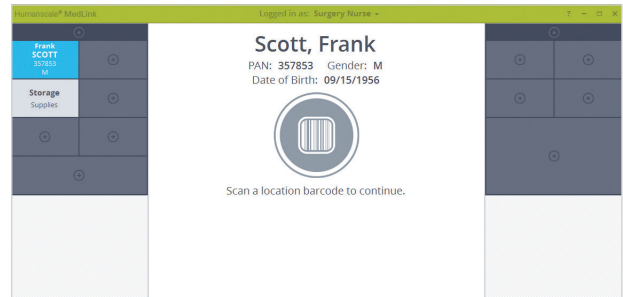


ELECTRONIC DRAWER LABELING

- Drawers can be assigned to an individual patient through electronic labeling
- HL7 integration allows patient labels to be selected from the existing hospital patient list rather than manual entry
- Drawers may be labeled as storage drawers when transporting non-patient specific supplies or consumables
- Search for patients by department
- RFID-embedded technology assists in locating each patient drawer — to the specific cart, cassette and caregiver

BAR CODE SCANNING

- Assigned drawers can be opened automatically by scanning a patient wristband barcode
- Drawers can be opened and closed individually
- Barcode Scan Locations (at the ADC or in a medication store room) ensure drawers can only be accessed in a designated area when filling medication
- Override events (opening drawers without scanning the patient's wristband or location barcode) are captured via reporting software

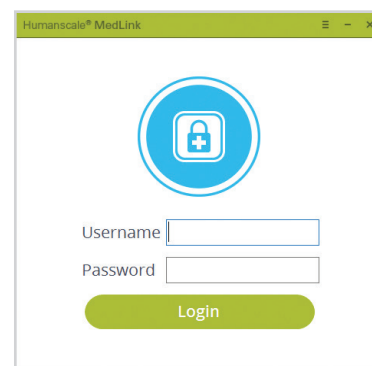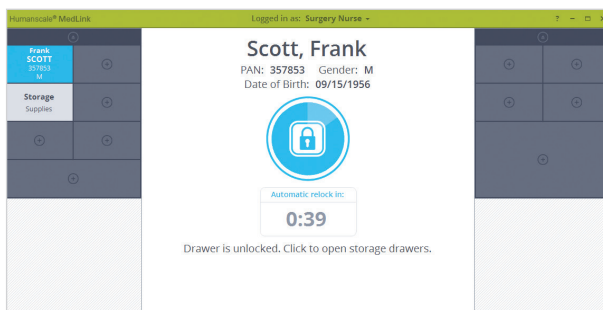

AUTO-CLOSE

- If patient drawers are left open, the system will automatically close and lock the drawer, preventing unauthorized access to medication if the caregiver is called away

IT INTEGRATION

- Microsoft Active Directory integration, with option to set-up single sign-on capabilities
- Access provided through secure username and password system


Effortless Mobility and Height Adjustment

MedLink incorporates features of Humanscale's TouchPoint™ Mobile Technology Carts.

These features provide MedLink users with unmatched mobility, exceptionally easy and comfortable use, lasting durability and high performance.


POWER TRACK™

With Power Track™ steering, caregivers can effectively navigate tight spaces, corridors and turns with one-hand operation for care on the go. Power Track is conveniently actuated with the effortless squeeze of a button on the underside of the cart handle. A fifth, center wheel acts as an axis point, allowing the cart to make sharp turns with one hand.

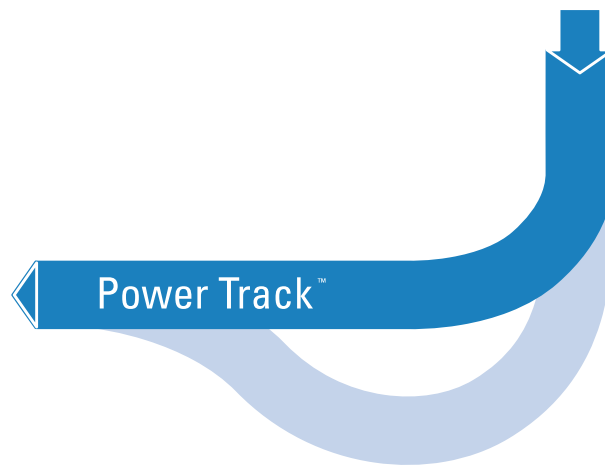

AUTO FIT™

Revolutionary Auto Fit™ technology instantly adjusts to each caregiver's height with a single touch and supports caregivers' increased reliance on computers in very short intervals.

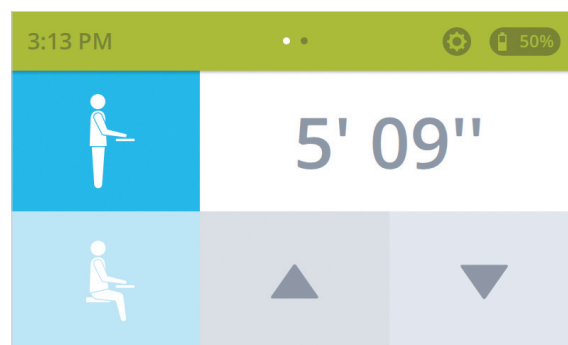

Which MedLink is Right for You?

Feature	MedLink Lite	MedLink Pro
Type of healthcare environment	Long term care environment with low patient turnover	Acute care environment with high patient turnover
Typical cart usage	A cart is assigned to a single nurse for the duration of their shift	Carts are shared by several nurses over the course of the day
Security requirements	A single PIN code for accessing medication is acceptable	Access to medication needs to be controlled on a per user basis
Reporting requirements	Detailed records of cart usage and access to medication are not needed, or is captured by other systems and processes	Detailed reporting and tracability of cart usage and access to medication is required

FEATURE COMPARISON

Feature	MedLink Lite	MedLink Pro
Electronically open individual drawers	✓	✓
Electronically open all drawers at once	✓	✓
Drawer Auto-Close (prevents drawers being accidentally left unlocked)	✓	✓
4-digit PIN code controls access to drawers	✓	✗
Supports multiple users/passwords	✗	✓
User-specific password controls access to drawers	✗	✓
Centrally manage Users, Groups, Departments	✗	✓
Reporting – tracking each step throughout workflow	✗	✓
Advanced login/single sign on options (requires Microsoft Active Directory integration)	✗	✓
Patient ID (requires HL7 integration)	✗	✓



1300 HPAUST | info@hpaust.com | www.hpaust.com

Humanscale[®]
HEALTHCARE

1-800-400-0625 | www.HumanscaleHealthcare.com

www.productportal.humanscale.com/medlink