### TADLI The Australian Dermal & Laser Institute


@ MEICET "


# **OMEICET**

# THE FUTURE OF

The Meicet MC10 Skin Scanner utilises advanced image analysis and processing technology to provide comprehensive skin diagnostics.

This innovative device takes images under five spectral photography modes. The software then analyses these images using algorithmic techniques to generate a total of 12 detailed images, revealing sensitivity, pigmentation, pores, wrinkles, hyperpigmentation and porphyrins.

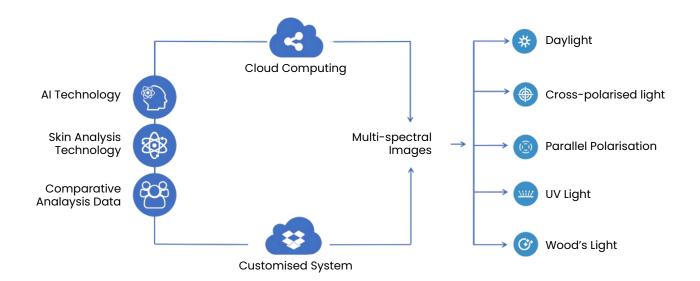

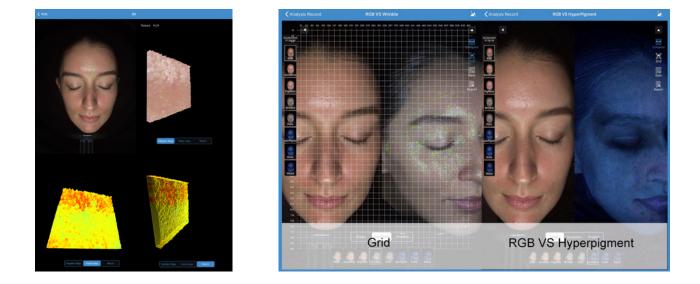

Distributed by The Australian Dermal & Laser Institute (TADLI)



## SPECTRAL PHOTOGRAPHY

Capture skin images under five lighting conditions for comprehensive analysis:

- 1. **RGB (Red, Green, Blue) Light** standard lighting to capture a high-resolution image of the skin
- 2. Cross-Polarised Light reduces surface glare and reflection, providing a clear view of deeper layers
- 3. **Parallel-Polarised Light** illuminates the skin at specific angles to enhance visibility of skin texture
- 4. UV Light reveals sun damage that are not visible under normal lighting
- 5. Woods Light a specific type of UV light to detect dermatological conditions















1.RGB Photo 2.Sensitivity & Heat Map 3.Spots & Spots Ageing 4.Wrinkles 5.Pores 6.UV Spots & Monochrome 7.Woods

Distributed by The Australian Dermal & Laser Institute (TADLI)

www.dermalandlaser.com.au

## MC10 KEY FEATURES

### **User-Friendly Design**

Equipped with a chin rest, magnetic blackout cover and a removable iPad bracket.

### **Comprehensive Reports**

The software generates a report including skin age, ratings of sensitivity, pigmentation, wrinkles, pores, UV spots and porphyrins, and detailed analyses of possible causes and management for each condition.


### **Product Recommendations**

Upload retail products to the software, allowing the scanner to recommend tailored solutions based on each client's unique skin profile. From the recommended products, select which products to add to the report.

### **Client Communication**

Easily email a report with images, score results, notes and recommended products directly to your clients, all through the app.

Test Decest				
-	Test Report			Test Report
	🥹 📀			Larimedical Melanze Cream 50ml 170.44 Symptom UV Spots(Mid, Moderate, Severe) Tips Use AM daily, after cleansing and serums.
*				Larimedical Hydrapeel Serum 15ml 106.62 Symptom Sensibirity(Mid, Moderate, Severe), Spots(Mid, Moderate, S. Tips Use PM daily, after cleansing, finish with cream.
				Larimedical Comedopeel Serum 15ml 133.16
				Symptom Spots(Mild, Moderate, Severe), Porphyrins(Mild, Moderate, S. Tips Use PM daily, after cleansing, finish with cream.
esuits Comprehensive scores: 38 Detection Herms		Bin Age: 64		
Detection Items	Severity company	ed to peers		
Detection Items Sensitivity	Severity compare 99.00%	nd to peers * * * \Severe		
Detection Items Sensitivity Spots	Severity compare 99.00% 49.53%	d to peers ★ ★ ★ \Severe ★ ★ ★ \Moderate		
Detection Items Sensitivity Spots Wrinkles Pores UV Spots	Severity compare 99.00% 49.53% 55.00% 61.73% 62.02%	d to peens ★ ★ ★ \Sovere ★ ★ ★ \Moderate ★ ★ ★ \Moderate ★ ★ ★ \Moderate ★ ★ ★ \Moderate		
Detection Hems Sensitivity Spots Whikles Pores	Severity compare 99.00% 49.53% 59.00% 61.73%	id to peens ★ ★ ★ \Severe ★ ★ ☆ \Moderate ★ ★ ☆ \Moderate ★ ★ ☆ \Moderate		

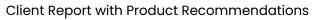

Raising the standard of skin education & technology and bridging the gap between beauty, dermal & health care practitioners, ensuring an integrated approach.

Award winning, bespoke

device training courses

& clinical coaching.

IPL & Laser

LED & LLLT

Skin Needling Electrolysis Diathermy

# ABOUT TADLI

Established in 2007, The Australian Dermal & Laser Institute (TADLI) offer diverse and comprehensive aesthetic device training for various manufacturers and direct to clinics. We believe that "best practice" is not negotiable and safe, superior results are the minimum standard. This is key to having highly credible staff members, ensuring greater client satisfaction, and guaranteeing a successful business.

TADLI are proud to offer an array of aesthetic devices that are both cost effective and results driven. Our strong focus on education and maximising treatment results, makes us a trusted partner to ensure clinics receive comprehensive training and support.

### Why Partner with TADLI?

### Award-Winning Training:

Every device includes award-winning training, providing users the expertise to maximise potential.



### Continuous Clinical Support:

We stand by our partners with ongoing clinical support, providing insights, solutions and advice whenever required.



### Trust and Integrity:

Our core values emphasise building enduring relationships based on trust, making us a dependable partner in your clinic's aesthetic journey.



### Promotion of Safe Practices:

With our extensive training and continuous support, clinics can ensure they're always delivering treatments that are safe and compliant.



### Authentic Female Leadership:

Founded and run by women who've thrived in their own aesthetic clinics, we prioritise genuine care and industry expertise over mere profit-making.

# MC10 KEY SPECIFICATIONS

Spectra	RGB Light, Cross-polarised Light, UV Light, Parallel-polarised Light, Wood's Light
Lighting Technology	Solid-State LED
Adapter	Input: AC100-240V, 50-60Hz, 1.5 max Output: DC 24V, 2A
Power Consumption	40 W
Net Weight	8 kg
Gross Weight	10 kg
Main Unit Dimensions	400x430x550mm
Package Dimensions	530x500x570mm
iPad Compatibility	5th-10th generation iPads





Distributed by The Australian Dermal & Laser Institute (TADLI) info@dermalandlaser.com.au | (03) 8001 5004 | www.dermalandlaser.com.au