Embedding Dispensing Console

	-		
Voltage/Frequency		220 V ± 10%, 50 HZ / 110 V ± 109	6,60 HZ
Nominal power		1200 W	
Working temperature		from ambient to 90℃ Increments 1	$^{\circ}$
Capacity			
Paraffin reservoir		6L	
Cassettes/molds tray		1.4L (126 cassettes)	
Dimensions			
Width		560mm	
Depth		640mm	
Height		400mm	
Weight		about 26 KG	

Cryo Console

Voltage/Frequency	220 V ± 10%, 50 HZ / 110 V ± 10%, 60 HZ
Nominal power	330 W
Working temperature	15℃~-15℃ Increments 1℃
Dimensions	
Width	390mm
Depth	560mm
Height	390mm
Weight	about 24 KG
Capacity	70 cassettes

Thermal Console

Voltage/Frequency	220 V ± 10%, 50 HZ / 110 V ± 10%, 60 HZ
Nominal power	335 W
Working temperature	from ambient to 75℃ Increments 1℃
Capacity	
Cassettes/molds tray	2.5L (252 cassettes)
Dimensions	
Width	300mm
Depth	540mm
Height	390mm
Weight	about 7.5 KG





AMOS SCIENTIFIC PTY.LTD.

Address: 7 Gabrielle Court, Bayswater North VIC 3153 ,Melbourne, Australia Tel: +61-410190188 Fax: +61-3-98081503 Email: info@amos-scientific.com





Embedding Dispensing Cosole

- Large paraffin tank capacity: 6 Liter
- The lamp brightness can be easily adjusted, according to ambient light brightness
- The large LCD screen is based on user– friendly design
- Automated "Wake up" function can be easily programmed

Cryo Console

- The paraffin block can be cold down rapidly by the stable and environmental friendly compressor
- The temperature range is between −15 to 15° C
- 70 cassettes can be loaded on this large work surface at the same time
- User-friendly operation interface



Characteristics



Working surface

- The extra large working surface allows the pathologist to process more samples simultaneously
- The paraffin flow can be precisely controlled by the finger touch pad or the foot pedal
- The paraffin flow rate can be controlled by the knob
- The "built-in" forceps heating station can heat 6 forceps simultaneously

FOOT PADEL



Thermal Console


- The Thermal console is compatible with all types of cassettes, embedding molds and processed tissues specimens
- The large capacity: 2.5L (252 cassettes)
- The temperature range can be set from ambient to 75°C
- This thermal console can also be used as a water bath





- The full information LCD screen shows current/setting temperature of paraffin tank, left/right tray, platform surface and forceps heating station respectively
- The peltier cold spot's temperature range is from ambient to −5°C



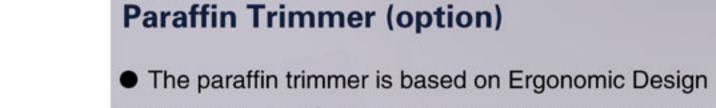

- The heated surface allows the users to remove excess paraffin easily
- The temperature range on surface can be set between ambient to 90° C
- The user friendly interface allows users to operate it easily and safely

