


Westlab Ductless Fume Cupboards

Codes 664-970, 664-971

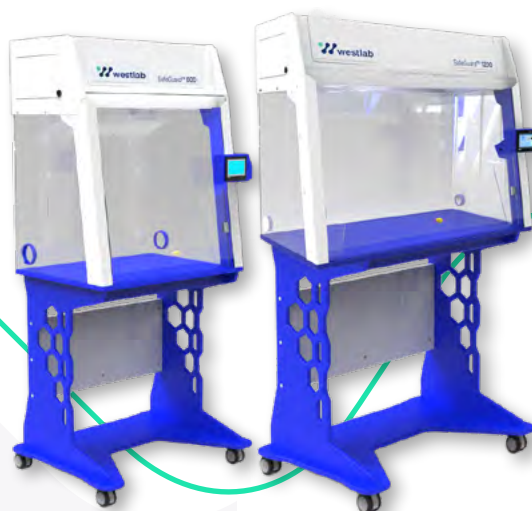

Westlab Ductless Fume Cupboards

Westlab's Ductless Fume Cupboard is a portable and economical unit designed to protect the operator from toxic gases, fumes, and noxious particulates.

With a compact and modern design, the filtered fume cabinet features all filters and electrics above the workstation to eliminate the risk of damage due to chemical spillage. Additionally, accidental spills are contained within the liquid tight base, without deducting from the workspace.

To ensure proper maintenance, a low airflow, gas, and filter replacement alarm signals the user when necessary. Furthermore, the ductless fume cabinet exceeds the requirements of AS/NZS 2243.9:2009 Safety in Laboratories – Recirculating Fume Cabinets.

Designed and manufactured in Australia, the fume hood requires no ductwork. Plugging directly into a standard 240V 10AMP electrical outlet, the ductless fume hood offers minimal installation.

Ideal for classroom environments and smaller laboratories.

Product Features

- Ductless Design, Sliding Sash
- Deep Base for Isolation of Accidental Spills
- Electrics & Filters Located Above Workstation
- Low Airflow, Gas & Filter Replacement Alarm
- Carbon Filtration Suitable for Chemicals
- Lighter Than Traditional Metal Units
- Internal Baffle for Uniform Airflow
- Manufactured In Australia
- Operator Safety Controller
- Corrosion Proof Material
- LED Lighting, Cord Ports
- Minimal Installation

PRODUCT CODE

664-970
664-971

Australia

1800 358 101
sales.au@westlab.com
www.westlab.com.au

Canada

1 877 822 1455
sales.ca@westlab.com
www.westlab.com

New Zealand

+64 9 553 3677
sales.nz@westlab.com
www.westlab.co.nz

United States

646 568 5391
sales.us@westlab.com
www.westlab.com

Product Specifications

Sliding sash

Code	664-970	664-971
Speed Selection	Variable speed	Variable speed
Low Airflow Alarm	Audible/visible	Audible/visible
Volume of Filtered Air	80 to 100L per second	120 to 140L per second
Linear Face Velocity	0.50M per second	0.50M per second
Top & Base Material	Polypropylene	Polypropylene
Front Panel Material	PETG	PETG
Side & Back Material	PETG	PETG
Base/Worktop Material	HDPE	HDPE
Voltage	110 to 240V	110 to 240V
Current	0.5A at 240V	1.0A at 240V
Power Consumption	120W	240W
Internal Dimensions	790mm W x 685mm H x 500mm D	1194mm W x 685mm H x 500mm D
External Dimensions	800mm W x 990mm H x 600mm D	1200mm W x 990mm H x 600mm D
Weight [Filters Included]	56.6kg	79kg

Accessories/Attachments

Code	664-971
Description	Ductless Fume Cupboard 800mm Trolley
Dimensions	800mm W x 600mm D x 880mm H

Code	664-973
Description	Ductless Fume Cupboard 1200mm Trolley
Dimensions	1200mm W x 600mm D x 880mm H

Code	664-976
Description	General Purpose Filter 600mm
Dimensions	610mm L x 305mm W x 76.2mm H
Weight	9kg

Code	664-977
Description	General Purpose Filter 450mm
Dimensions	457mm L x 305mm W x 76.2mm H
Weight	5.6kg

Australia
1800 358 101
sales.au@westlab.com
www.westlab.com.au

Canada
1 877 822 1455
sales.ca@westlab.com
www.westlab.com

New Zealand
+64 9 553 3677
sales.nz@westlab.com
www.westlab.co.nz

United States
646 568 5391
sales.us@westlab.com
www.westlab.com

Ordering Information

Code	Description
664-970	Ductless Fume Cupboard 800mm
664-971	Ductless Fume Cupboard 800mm Trolley
664-972	Ductless Fume Cupboard 1200mm
664-973	Ductless Fume Cupboard 1200mm Trolley
664-976	General Purpose Filter 600mm
664-977	General Purpose Filter 450mm

Australia

1800 358 101
sales.au@westlab.com
www.westlab.com.au

Canada

1 877 822 1455
sales.ca@westlab.com
www.westlab.com

New Zealand

+64 9 553 3677
sales.nz@westlab.com
www.westlab.co.nz

United States

646 568 5391
sales.us@westlab.com
www.westlab.com