

DERM 101/102 HIGH FREQUENCY DESICCATOR

The Derm 101 and Derm 102 are economical, 10 watt high frequency desiccators, thoughtfully designed to give practitioners the ability to perform in-office, minor skin procedures such as the removal of skin tags, benign lesions and premalignant lesions. The Derm 101 High Frequency Desiccator offers ten watts of monopolar power, adjustable in 1/10th watt increments, to desiccate and fulgurate tissue. The power can be adjusted using the 3-Button Power Control Handpiece (A902) or using the rotary dial on the unit. The unit features a large, bright, LED display for excellent visibility across the procedure room. For the practitioner utilizing bipolar electrosurgery, Bovie offers the Derm 102 High Frequency Desiccator with 10 watts of bipolar power, adjustable in 1/10th watt increments.

Style # DERM 101 & 102

Dimensions

Height: 208mm (8.2 in.)
 Width: 184mm (7.2 in.)
 Depth: 69mm (2.7 in.)
 Weight: < 1.4 kg (< 3 lbs.)

Output Power

Up to 10 watts (0-10 watts in 1/10th watt increments on both monopolar and bipolar)

Main Voltage

Universal Power Supply
 Derm 101/102: 100 to 240 VAC±10%

Main Frequency: 50 – 60 Hz

Main Current: 320mA Max.


Warranty: **FOUR YEARS**




DERM 101 on A813 Table-top Stand*



DERM 102 on A813 Table-top Stand*

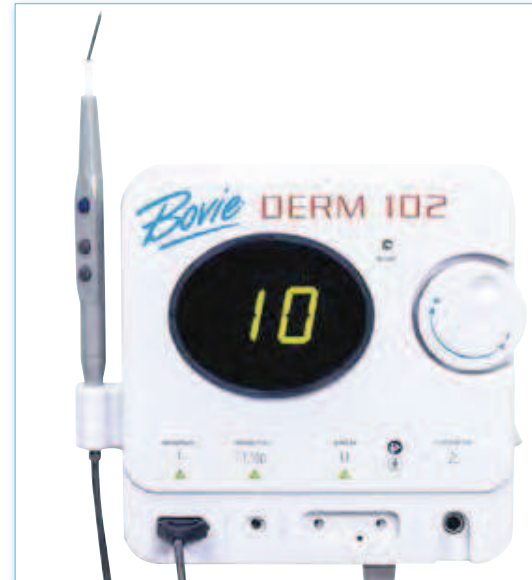

Power Cord

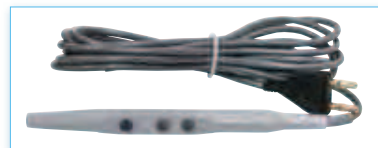

A902



A804, A805



A806, A807



A910



A837

* Sold Separately

Standard Equipment

Derm 101 or Derm 102 unit

110 VAC medical-grade power cord (09-064-001)

1 three-button power control autoclavable handpiece (A902)

2 non-sterile, disposable handpiece sheaths (A910)

5 non-sterile single-use sharp dermal tips (A804)

2 sterile single-use sharp dermal tips (A805)

5 non-sterile single-use blunt dermal tips (A806)

2 non-sterile single-use Derm-Elite® blunt dermal tips (A806DE)

2 sterile single-use Derm-Elite® blunt dermal tips (A807DE)

1 wall mount kit (A837)

1 user's guide CD (MC-55-229-001)