



TIDISHield®

Now with Point-of-Use Dispensing Options



*TIDISHield® Assemble 'n Go™ Face Shields
PROTECTION and COMFORT*

PROTECTION IS ESSENTIAL – TIDISHIELD® ASSEMBLE 'N GO™ FACE SHIELDS ARE THE ANSWER

Healthcare professionals work under hazardous and ever-changing conditions. What healthcare workers need is superior face protection from the splash and splatter of blood-borne pathogens and other potentially infectious materials. TIDISHield® Assemble 'n Go™ Face Shields are the answer.

TIDISHield® Assemble 'n Go™ Face Shields are a cost-effective face protection product. Face shields are lightweight, cool and offer comfortable protection that can be worn over prescription eyewear and facemasks.

THE NEED FOR PROTECTIVE FACEWEAR

- The mucous membranes of the mouth, nose, and eyes are susceptible portals of entry for infectious agents, as can be other skin surfaces if skin integrity is compromised.¹
- The wearing of masks, eye protection, and face shields in specified circumstances when blood or body fluid exposures are likely to occur is mandated by the OSHA Bloodborne Pathogens Standard.¹

COMFORT IS A PRIORITY

- Lightweight design can be worn for long periods of time (< 1.5 oz)
- Adjustable face shield provides flexibility for the worker to move into optimal position
- Includes a nose piece for comfortable and secure fit, even with glasses
- Maximizes airflow to keep the face cooler as compared to blue foam face shields
- Good for sensitive skin and avoiding skin breakouts
- Patented Design

TIDI - Assemble 'n Go Face Shield Dispensers

Provides protective eyewear when and where you need it.

Designed to maximize airflow and prevent fogging- weighs less than 1.5 ounces! Provides over-the-top, side and front Face protection. Disposable, optical grade, anti-static replacement shields. Includes a clear, snap-in nosepiece for support.



2209

TIDI FACE SHIELDS DISPENSER SET

- * 20 Visors,
- * 20 Nose Pieces
- * 40 Shields per Dispenser



2211-80DSP

TIDI FACE SHIELDS DISPENSER SET

- * 20 Visors,
- * 20 Nose Pieces
- * 80 Shields per Dispenser
- * Assorted Colors



2210-100

TIDI FACE SHIELDS

- * 100 Shields
- * Use with Visors 13" x 7"



2210-40

TIDI FACE SHIELDS

- * 40 Shields
- * Use with Visors 13" x 7"

TIDI VISORS

- * (2211-100)
- * Assorted Colors
- * 100 (10 Bags of 10 Visors)



7005

TIDI CLEAR FACE SHIELD DISPENSER HOLDERS

- * (7005)
- * 2 Holders per Case
- * 8.5"W x 4.75"D x 8.125"H


ASSEMBLY INSTRUCTIONS

For ultimate under face protection, place visor on the bottom hole, for more room, place visor on the top shield hole.

(1) Insert nosepiece tab into visor slot, so nosepiece slant away from front of visor. (Nosepiece may be removed to wear over eyeglasses)



(2) Starting on one side of the visor, first mount the face shields onto the visor pin (1), then insert the visor's support knuckle (2) through the face shield slot

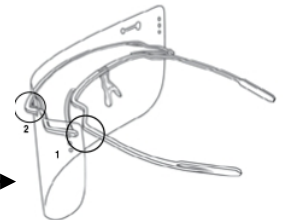

(3) Repeat step 2 on the opposite side of the visor, pulling the face shield tight and securing it to the opposite support knuckle and pin.



Selection of mounting holes provides 3 options for adjusting the face shield angle. Top hole will position the face shield farther from chin, and bottom hole will position it closer to chin.



Majac Medical Products Pty Ltd

Ph: +617 3265 6355
 Fax: +617 3865 2729
 sales@majacmedical.com.au
 www.majacmedical.com.au



FS 540017

Also Available From: