

SCARLET

AUTHENTIC

Minimal invasive fractional micro-needle
RF with micro-needle electrodes system.

The only skin needling RF device with
the patent "Na Effect" technology.

SKIN TIGHTENING & LIFTING

ACNE SCAR REMOVAL

PORE & STRETCH MARK REDUCTION

somaTM
Lasers Australia

FDA TGA ISO CE



Scarlet S

The first true fractional micro-needle RF

Scarlet S offers total skin rejuvenation treatments. It is the safest treatment in the field of RF systems, with less downtime and side effects. Clients can resume daily activities after treatment and there is no bleeding or sensitivity as experienced with other micro-needle treatments. It targets virtually any skin issue. Treatments can be done on the face and body.

It provides direct heat to the epidermis and dermis by precise depth control from 0.5mm to 3.5mm. Its efficient fractional radio frequency offers fractional and continuous treatment patterns.

Scarlet S micro-needling RF stimulates collagen production, which gives skin a healthy, radiant glow. RF also works to tighten the skin, which can eliminate fine lines in certain areas. RF micro-needling can also reduce the appearance of some scars, such as acne scars.



“Na Effect”

What is the patented “Na Effect”?

Dr. Jongju Na, from South Korea, pioneered the “Na Effect” that delivers energy using microneedling in a unique way. This patented solution safely and effectively tightens skin in fewer treatments. The “Na Effect” creates a wound like response in one pass – competing devices take up to three passes at variable depths to achieve similar results. This advantage, combined with superior speed and finely tuned presets, makes Scarlet S RF treatments more comfortable and fast.



The energy follows the path of least resistance and the tip is where the path of least resistance lies. This flow causes the energy to dissipate from there and work its way to the top, this distribution of energy creates a unique teardrop effect. Scarlet S RF is intended for use in dermatologic and general surgical procedures for electrocoagulation and homeostasis.



Pictured: Dr. Jongju Na


Difference of Non-Insulated Vs Insulated

Non-insulated Scarlet S micro-needle delivers RF energy fully and safely on epidermis and dermis resulting in no bleeding during the treatment on the patient's skin.

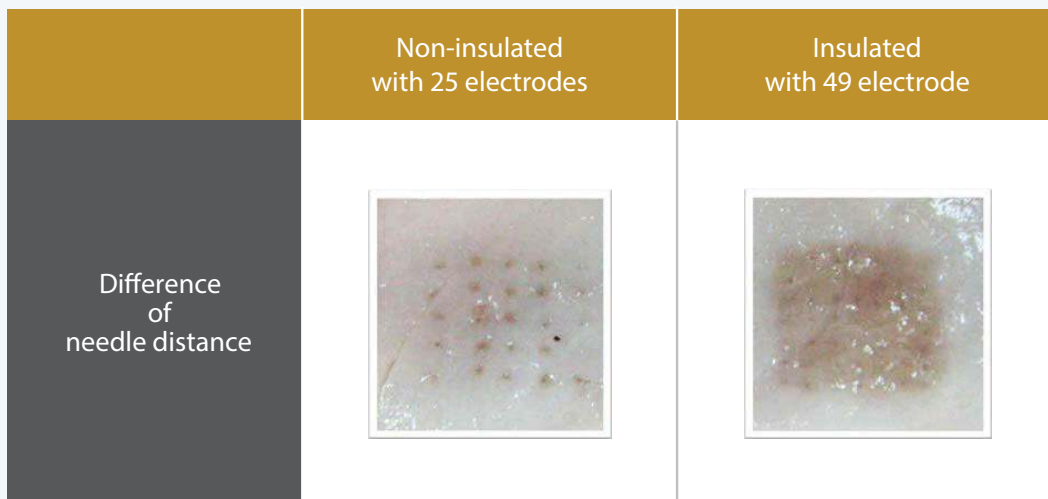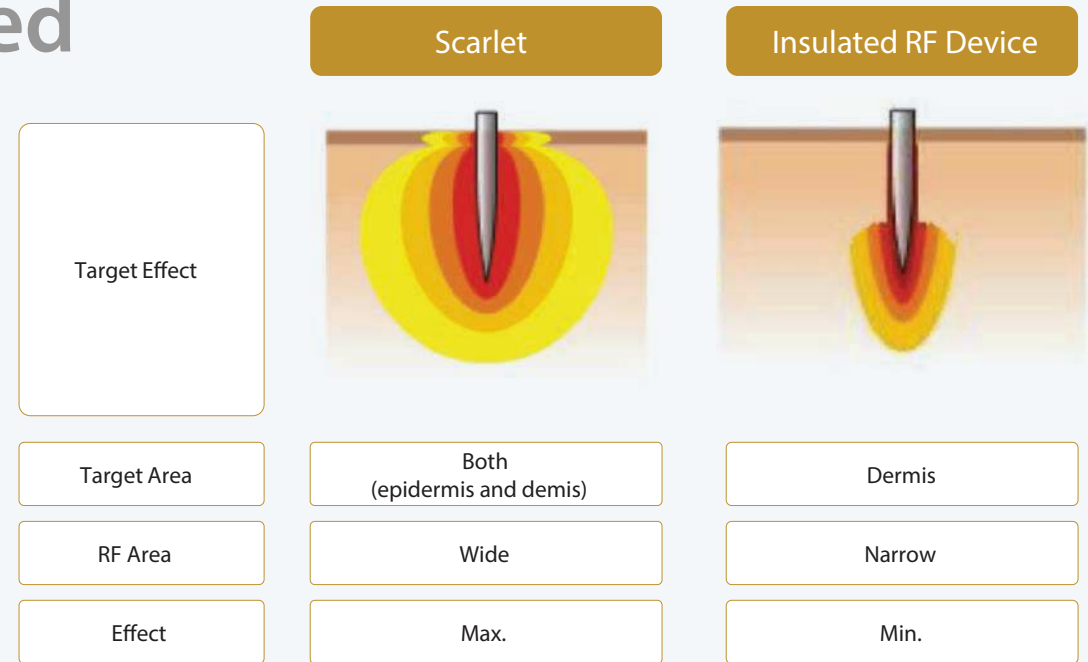


Enough and high RF energy can be applied to the target areas without any risk of energy interference caused by overlapping of RF energy in close distance of treatment. There is no fear of skin burn erythema and edema.

Treatment

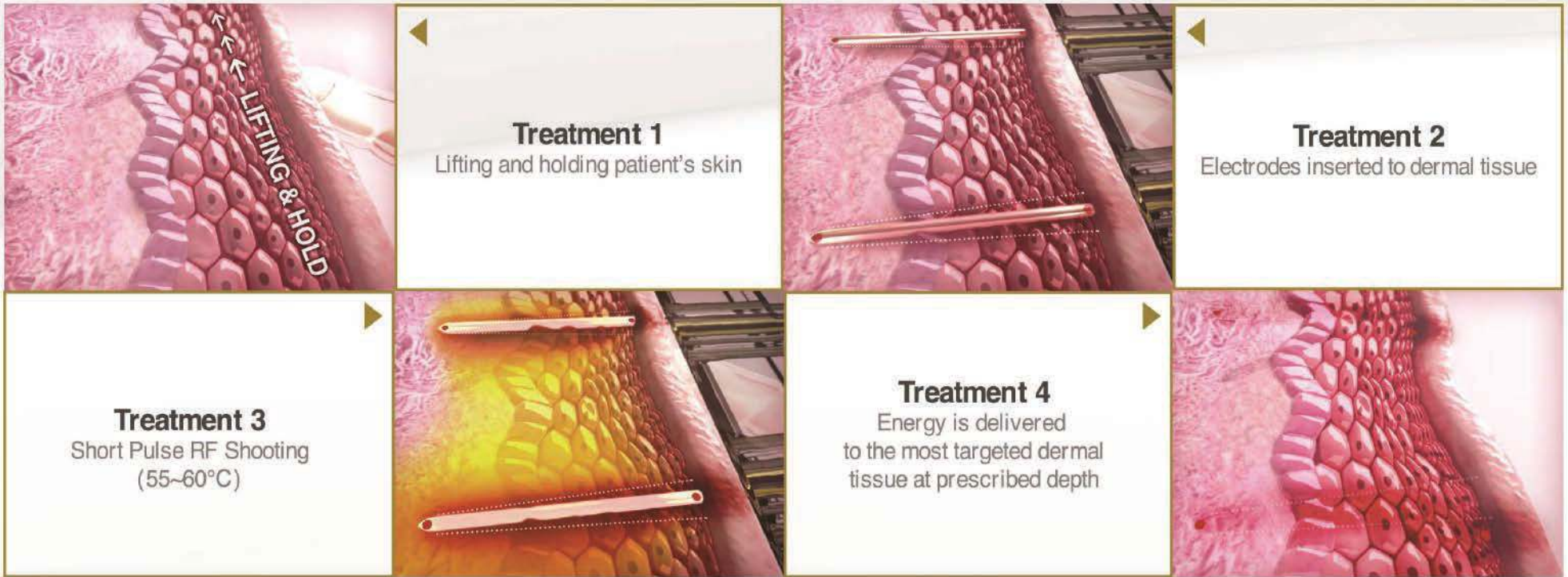

Result

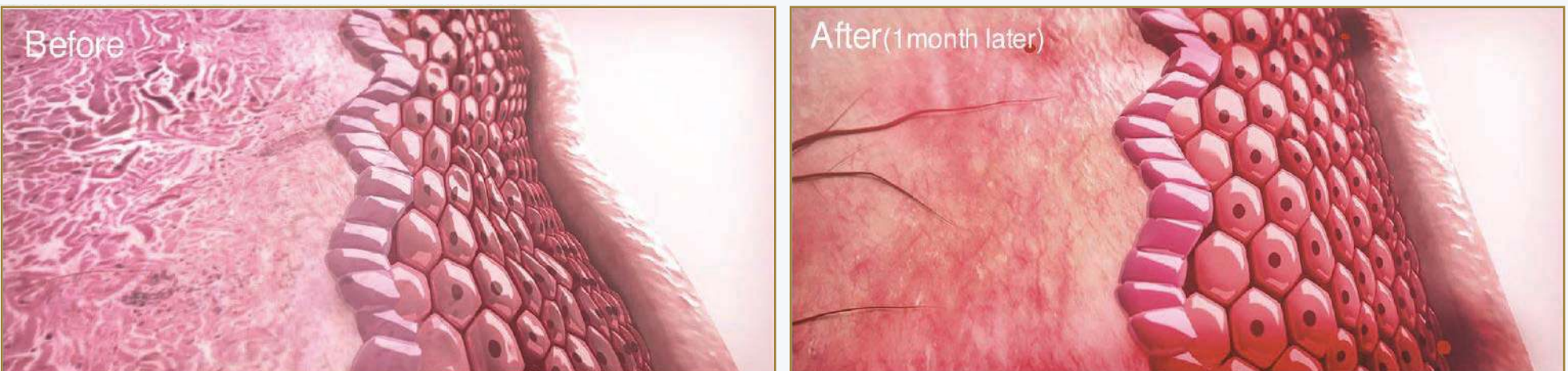

CLINICAL DATA



“Quality in a product is not only what you put into it.
It is what the customer gets out of it.”

~ Soma Lasers Australia



Unmatched, incomparable device made with the therapist in mind

Soma devices are made with the highest quality engineering working together with Germany, Korea, Japan the USA and proudly designed in Australia for fantastic durability.

Warranty

We back up our promise at Soma Lasers Australia, with a 2 year warranty on all our new medical aesthetic devices.

Scarlet S Specifications

System type	Bipolar RF (Radio Frequency)
Frequency	2MHz
Total RF treatment level	Min. (Level 1) to Max. (Level 10)
Tip depth	Micro-needle electrodes, 0.5-3.5mm (adjustable)
Spot size	10x10mm (25 pin)
Display	Dual LCD (colour touch LCD)
Dimensions	550 x 350 x 1220 (mm)
Voltage	AC 100-260V, 50/60Hz
Weight	30kg

* All specifications are subject to change without notice

* Please read precaution and instructions before use

Equip your salon
or clinic with



Our demonstration suite is situated at:
3/10 Atherton Road, Oakleigh Vic 3166

M: +61 403 744 648
E: support@somalasersaus.com.au
W: www.somalasersaus.com.au