

ClinicalPRO™
SKINCARE + EQUIPMENT + EDUCATION



755nm 808nm 1064nm Three Wavelengths Diode Laser

TRINITY®

Liftoning





TRINITY[®] Liftoning

TRINITY is a medical device that emits lasers of three different wavelengths simultaneously. LIFTONIG is the representative treatment of Trinity which has skin lifting, tightening and brightening with wide range of penetration for better photo thermal effects and non-invasive skin procedures effects through photothermal accumulation.

WHY

TRINITY[®]

Liftoning



Three Wavelengths Technology



Three Contact Cooling System



Three Treatment Modes



Three Treatment Results
Lifting / Tightening / Brightening

THREE CONTACT COOLING LEVEL


Handpiece of Trinity has three levels of contact cooling for different patient's skin conditions to support comfortable treatments with large irradiation spot size (12×20mm).

Sapphire Contact Cooling Tip



THREE TREATMENT MODES

Trinity offers three treatment modes for each patient's skin condition safely through the precise laser thermal energy.

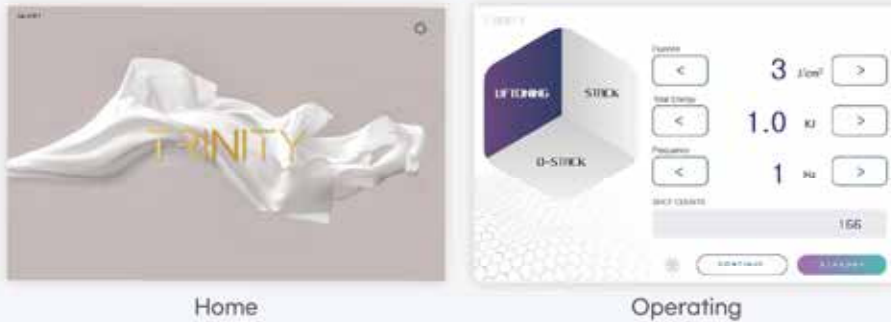


Three Lifting Modes

Various customized treatments with three lifting modes



THREE WAVELENGTHS TECHNOLOGY


Trinity is the first Korean three wavelengths diode laser operating on 755nm, 808nm, 1064nm simultaneously. Trinity's unique technology allows it to effectively target various depths within the skin through its three wavelengths technology

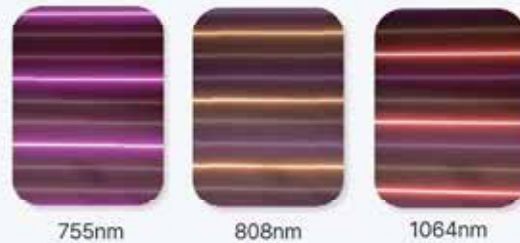


755nm Alex	Simultaneous
808nm Diode	Irradiation of
1064nm Nd:YAG	Three wavelengths

Versatile Laser

When the energy of a multiple wavelengths laser with different skin absorption depths reaches human tissue, the chromophore of the living tissue absorbs the energy and exhibits a selective broad-area effect, thereby showing effectiveness in pigmented lesions and the heat applied to the irradiated area allows to cut, destroy and remove tissue.

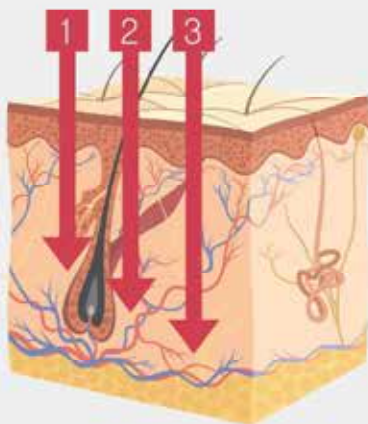

1 755nm 2 808nm 3 1064nm

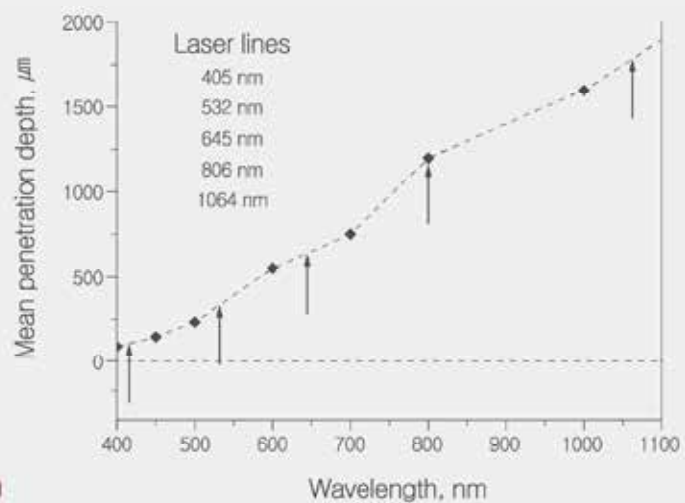








It is now the era of
TRINITY[®] Liftoning

Feel The Confidence With Three Wavelengths



Benefits of TRINITY[®] Liftoning

-  10Hz balanced laser output
-  No anesthesia
-  Contact cooling system
-  Non-invasive procedure
-  10 minutes treatment
-  12X20mm irradiation area
-  Simultaneous irradiation of three wavelengths

SPECIFICATIONS

Laser Type	Diode Laser
Wavelength	755nm+808nm+1064nm
Spot Size	12X20mm
Power	800W
Fluence	Max. 60J/cm ²
Pulse Duration	Max. 300ms
Pulse Repetition Rate	1-10Hz
Display	12.1" LCD
Auto Management	WiFi
Basic Part	12X20mm Handpiece
Electrical Requirement	220~230VAC, 15A, 50Hz/60Hz
Dimensions	380mm X 523mm X 346mm / 40kg



TRINITY[®]
Liftonging



ClinicalPRO[™]
SKINCARE + EDUCATION + EQUIPMENT

www.clinicalpro.com.au

For enquiries:
1800 628 999
ask@clinicalpro.com.au

ClinicalPRO
2/10 Valente Close
Chermside QLD 4032