

Floating X-Ray Tube / Table

Features

Our floating X-Ray table now enables you to accurately and simply get perfect positioning every time.


- The radiolucent x-ray table offers movement in every horizontal direction possible. After positioning, release the foot operated lock and the top locks noiselessly into position
- There is also now no need to disturb the patient when the cassette is changed or a grid is used, as both are held in trays underneath the table top
- As the x-ray unit is always automatically centered onto the cassette, beam alignment is assured every time





FLOATING X-RAY TUBE / TABLE

- Under surface self centering cassette tray
- Under surface grid tray
- Standard table top 1.6 m long
- Optional 2.0 m table top
- Optional adjustable counterbalanced x-ray unit height
- Optional Myelography plane
- Optional horizontal beam cassette holder
- Suitable for cassettes, CR imaging plates, and most DR panels
- Standard table top size (wxd): 1600 x 650 mm
- Optional table top size (wxd): 2000 x 650 mm
- Minimum room width for standard top: 2000 mm
- Minimum room width for optional top: 3000 mm
- Optional X-ray grid 45.1 x 43.8 cm. 6:1 103 LPI infinitive 17.76 x 17.25 " (pictured bottom right)

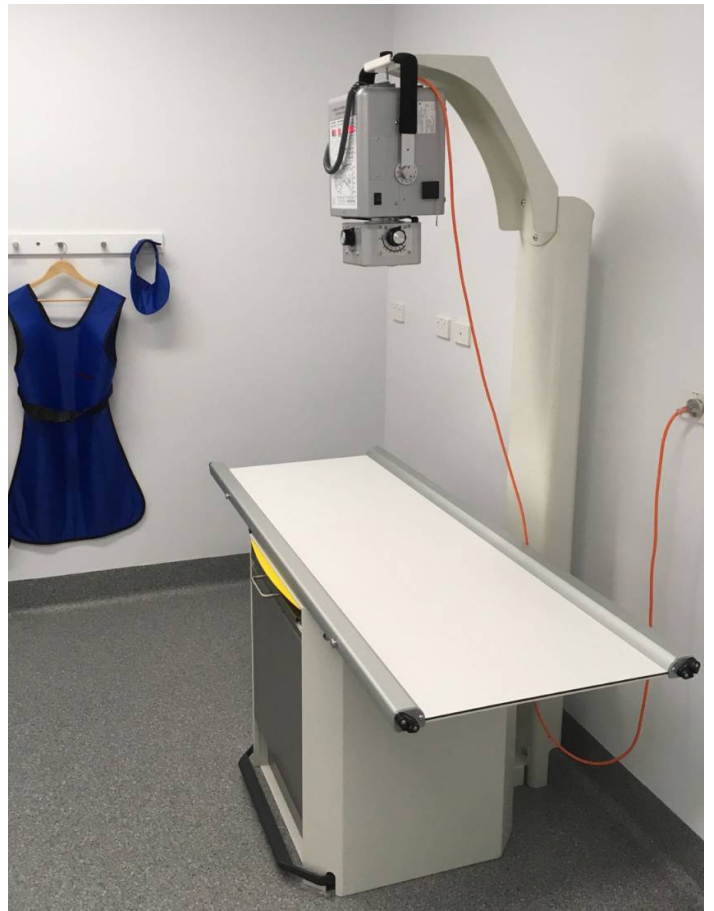

Atomscope X-Ray Unit not Included



330012.....X-Ray Table / Tube Stand with X-Ray Grid

330013.....X-Ray Table / Tube Stand without X-Ray Grid

For further information please call:



Southern office:

PO Box 4438,
17-19 Horne St,
Hoppers Crossing, VIC, 3029

Phone: 03 9360 9700
Fax: 03 9360 9994
Toll Free: 1300 785 405
Email: sales@dlc.com.au

Northern Office:

PO Box 833,
Unit 6/1 Lear-Jet Dr,
Caboolture, QLD, 4510

Fax: 07 5428 3777
Support: 1300 785 401
Website: www.dlc.com.au