

SPECIAL PRODUCTS FOR OEDEMA THERAPY

Compression bandaging with Juzo[®] SoftCompress



Juzo[®]. Freedom in Motion

Our customer service department – always there for you

Service with a capital S

To us, customer service means being reachable and always having a sympathetic ear. Our service teams will deal with your concerns with passion, flexibility and reliability. You can reach us by phone **MON to FRI** or 24/7 by email or fax.

Customer service also means expert advice. Take advantage of our more than 100 years of experience in compression therapy. Our friendly and dedicated staff will be happy to assist you with all your enquiries.

We look forward to your call!

Our services at a glance

- general information
- product advice
- order acceptance
- order status
- alterations and repairs

Customer Service



Edith Männer
Customer Service
Manager International

How to reach us

Juzo Customer Service International

Office hours 7:00 am - 6:00 pm

Phone +49 (0)8251 901 300

Fax +49 (0)8251 901 321

E-mail order@juzo.com

Juzo Customer Service UK

Office hours 9:00 am - 5:00 pm

Phone +44 (0)1382 82 66 20

Fax +44 (0)1382 82 66 41

E-mail order@juzo.co.uk

Juzo®

More than 100 years of experience
in compression therapy

Summary



Compression bandaging

Material and characteristics	04
How it works	05
Indications and contraindications	05



Bandage

Uni size	06
Standard size	06
Custom-made item	06



LymphAid and BandageAid

LymphAid	08
BandageAid	09



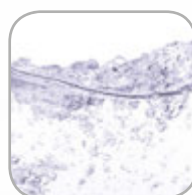
Compression roll and sheet

Compression roll	10
Compression sheet	10



Accessories

Breast pad and bandage	11
Genital pad	11
Knee and elbow joint pad	12
Malleo pad	13
Hand pad	13



Care instructions

Material composition	14
Washing and care instructions	14
Storage instructions	14

Compression bandaging

Material and characteristics



The vertically-stitched foam channels promote the lymph drainage.

Using Juzo SoftCompress products can significantly improve the treatment outcome. Compared to the often used foam or synthetic wadding, Juzo SoftCompress products are characterized by their particular skin friendliness and durability. All Juzo SoftCompress Bandages and Accessories are manufactured in Germany.

Characteristics

- especially skin-friendly surface made of cotton
- core made from foam
- vertically-stitched foam channels
- washable up to 60 °C
- breathable
- saves time when applying bandaging
- makes self-bandaging easier
- provides long-term cost savings



CE mark

Juzo SoftCompress products fulfill all the requirements of the Medical Products Act and European Council Directive 93/42/EEC.

How it works



Visible effect of the Juzo SoftCompress Bandage on a secondary leg lymphoedema.

The comfortable bandages can be worn under short-stretch bandages or compression stockings, and provide additional compression, cushioning and an ideal pressure distribution.

Juzo SoftCompress Bandages are employed as an effective treatment aid for oedema. The vertically-stitched foam channels create differences in pressure, thereby promoting lymph drainage, while at the same time preventing soft tissue fibrosis.

The bandaging aids avoid any constriction by the edges of the short stretch bandages and protect the bends of the joints. These characteristics and their effect are what makes Juzo SoftCompress so unique.

Indications

- lymphoedema
- phleboedema
- lipoedema
- fibrosis

Contraindications

- arterial insufficiency or degeneration
- deep vein thrombosis
- untreated congestive heart failure
- localized or systemic infection
- untreated malign cancer
- absent or severely impaired sensation

If there are any open wounds in the region being treated, these wound areas must first be covered with a suitable wound dressing or bandage before apply Juzo SoftCompress products.

Bandage

The comfortable Juzo SoftCompress Bandages can be worn under short-stretch bandages or medical compression stockings, and provide additional compression, cushioning and an ideal pressure distribution. The bandages are used for fast and effective bandaging and make self-bandaging easy for the patient. They are available in uni sizes, standard sizes and as custom-made items. Custom-made bandages are ideal for night time treatment.



Bandage lower leg in custom-made,
also available in standard sizes

6830



Bandage lower leg in uni size

6831



Standard and uni sizes can be cut
according to the cutting marks.



Bandage lower leg and thigh
in custom-made

6820

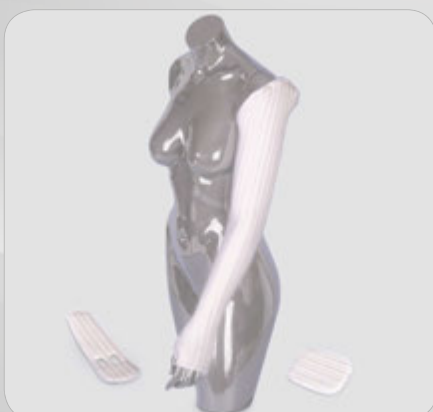


Bandage lower leg and thigh in uni size

6821

Uni size

Bandages in uni size have an open back side and can be fixed with elastic tapes or an elastic textile tube. Thereby the bandage can be adjusted individually on highly changing circumferences.



Bandage arm in custom-made,
also available in standard sizes

6822







Bandage arm in uni size



6823

For custom-made items please use the Juzo order forms for custom-made.

Juzo SoftCompress Bandages Leg

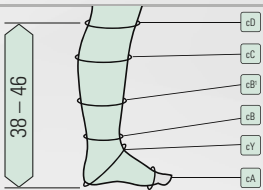
Description	Art. no.	Uni size	Standard size 1-5	Custom-made
 Bandage lower leg	6830		●	●
 Bandage lower leg incl. 2 malleo pads	6831	●		
 Bandage lower leg and thigh incl. 2 malleo pads 1 knee joint pad	6820			●
 Bandage lower leg and thigh incl. 2 malleo pads 1 knee joint pad 1 compression roll (10 cm x 4 m) 1 compression roll (12 cm x 4 m)	6821	●		

Juzo SoftCompress Bandages Arm

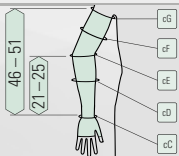
Description	Art. no.	Uni size	Standard size 1-6	Custom-made
 Bandage arm incl. 1 hand pad 1 elbow joint pad	6822		●	●
 Bandage arm incl. 1 hand pad 1 elbow joint pad	6823	●		

Size charts – standard sizes

Bandage lower leg • 6830

Lengths in cm	Circumferences in cm	1 x-small	2 small	3 medium	4 large	5 x-large
	cD	27 – 29.5	29.5 – 32.5	32.5 – 35.5	35.5 – 39	38.5 – 42
	cC	29 – 31.5	31.5 – 34.5	34.5 – 37.5	37.5 – 41	40.5 – 44.5
	cB'	23.5 – 25.5	25.5 – 28	27.5 – 30	30 – 33	32.5 – 35.5
	cB	18.5 – 20	20 – 22	22 – 24	24 – 26	26 – 28.5
	cY	28 – 30.5	29 – 31.5	31 – 34	33 – 36	35 – 38.5
	cA	18.5 – 20	20 – 22	22 – 24	24 – 26	26 – 28.5

Bandage arm • 6822

Lengths in cm	Circumferences in cm	1 x-small	2 small	3 medium	4 large	5 x-large	6 xx-large
	cG	30 – 33	32 – 35	34 – 37	36 – 39.5	38 – 41.5	40 – 44
	cF	28 – 30.5	30 – 33	32 – 35	34 – 37	36 – 39.5	38 – 41.5
	cE	25 – 27.5	27 – 29.5	29 – 31.5	30.5 – 33.5	32 – 35	33.5 – 36.5
	cD	23 – 25	24.5 – 27	26 – 28.5	27.5 – 30	29 – 31.5	30.5 – 33.5
	cC	14.5 – 16	16 – 17.5	17.5 – 19	19 – 20.5	20.5 – 22	22 – 24

LymphAid and BandageAid


Juzo SoftCompress LymphAid and BandageAid offer a complete solution for lymphological compression bandaging. Unlike conventional sets, they also include the innovative Juzo SoftCompress materials which can significantly improve the treatment outcome. LymphAid comes with a ready-to-use bandage which is also suitable for self bandaging. The more economical BandageAid contains compression rolls.

Juzo SoftCompress LymphAid



LymphAid leg small


6800 KL

Description	Art. no.	Set
 LymphAid leg small 1 bandage lower leg in uni size 20 gauze bandages 4 cm x 4 m (for compression of toes) 2 short-stretch bandages 6 cm x 5 m 2 short-stretch bandages 8 cm x 5 m 2 short-stretch bandages 10 cm x 5 m 4 short-stretch bandages 12 cm x 5 m 2 malleo pads	6800 KL	●



LymphAid leg large


6800 GR

Description	Art. no.	Set
 LymphAid leg large 1 bandage lower leg in uni size 1 compression roll 10 cm x 4 m 1 compression roll 12 cm x 4 m 20 gauze bandages 4 cm x 4 m (for compression of toes) 2 short-stretch bandages 6 cm x 5 m 2 short-stretch bandages 8 cm x 5 m 4 short-stretch bandages 10 cm x 5 m 6 short-stretch bandages 12 cm x 5 m 2 malleo pads 1 knee joint pad	6800 GR	●



LymphAid arm

6801

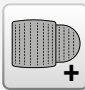
Description	Art. no.	Set
 LymphAid arm 1 bandage arm in uni size 20 gauze bandages 6 cm x 4 m (for compression of fingers) 2 short-stretch bandages 6 cm x 5 m 2 short-stretch bandages 10 cm x 5 m 3 short-stretch bandages 12 cm x 5 m 1 hand pad 1 elbow joint pad	6801	●

Juzo SoftCompress BandageAid



BandageAid leg small

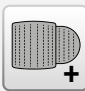
6810 KL

Description	Art. no.	Set
 BandageAid leg small 1 compression roll 10 cm x 4 m 20 gauze bandages 4 cm x 4 m (for compression of toes) 1 short-stretch bandage 6 cm x 5 m 1 short-stretch bandage 8 cm x 5 m 3 short-stretch bandages 10 cm x 5 m 2 malleo pads	6810 KL	●



BandageAid leg large


6810 GR

Description	Art. no.	Set
 BandageAid leg large 1 compression roll 10 cm x 4 m 2 compression rolls 12 cm x 4 m 20 gauze bandages 4 cm x 4 m (for compression of toes) 1 short-stretch bandage 6 cm x 5 m 1 short-stretch bandage 8 cm x 5 m 3 short-stretch bandages 10 cm x 5 m 4 short-stretch bandages 12 cm x 5 m 2 malleo pads 1 knee joint pad	6810 GR	●



BandageAid arm

6811

Description	Art. no.	Set
 BandageAid arm 1 compression roll 10 cm x 4 m 20 gauze bandages 6 cm x 4 m (for compression of fingers) 1 short-stretch bandage 6 cm x 5 m 1 short-stretch bandage 8 cm x 5 m 5 short-stretch bandages 10 cm x 5 m 1 hand pad 1 elbow joint pad	6811	●

The gauze bandages and short-stretch bandages in the LymphAid and BandageAid are quality products from the company Karl Otto Braun GmbH & Co. KG (KOB). They are also available separately.

Compression roll and sheet

Juzo SoftCompress Bandaging Materials are used for lymphological and phlebological compression bandaging. They are very supple on the skin due to the soft surface and also breathable.




Compression roll

6835

Juzo SoftCompress Compression Roll

The thin, very soft Juzo SoftCompress Compression Rolls already produce compression pressure if applied suitably. The compression rolls distribute the pressure of the additionally applied short-stretch bandages, thereby avoiding constrictions.

Description	Art. no.	Standard size
		1-4
 Compression roll size 1 – 5 mm, 8 cm x 4 m size 2 – 5 mm, 10 cm x 4 m size 3 – 5 mm, 12 cm x 4 m size 4 – 5 mm, 16 cm x 4 m	6835	●




Compression sheet

6836

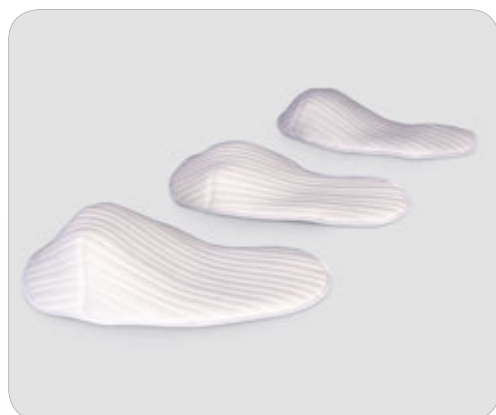
Juzo SoftCompress Compression Sheet

Juzo SoftCompress Compression Sheets are very thin and can be cut to size for use in various body regions.

Description	Art. no.	Standard size
		1-2
 Compression sheet size 1 – 65 x 65 cm size 2 – 125 x 100 cm	6836	●

Accessories

Juzo SoftCompress Pads can be worn under compression bandages. They are especially suitable for a targeted increase in pressure for particularly localized oedemas.




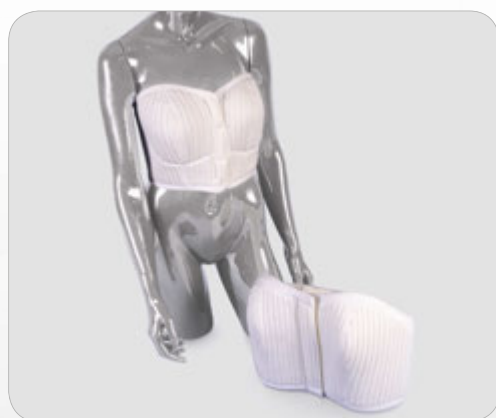
Breast pad

6837

Juzo SoftCompress Breast Pad

The Juzo SoftCompress Breast Pad serves as a comfortable pad for a compression bra or chest bandage, e.g. for breast or chest wall oedemas.

Description	Art. no.	Standard size	Custom-made
		1-3	
 Breast pad size 1 – cup A-B size 2 – cup C-D size 3 – cup E-F	6837	●	●




Breast bandage

6838

Juzo SoftCompress Breast Bandage

The Juzo SoftCompress Breast Bandage combines the positive characteristics of the Juzo SoftCompress material with those of a compression bra. It has an adjustable chest circumference. The breast bandage can, for example, be used for breast or chest wall oedemas.

Description	Art. no.	Standard size	Custom-made
		1-6	
 Breast bandage size 1 – cup A size 2 – cup B size 3 – cup C size 4 – cup D size 5 – cup E size 6 – cup F	6838	●	




Genital pad women

6839

Juzo SoftCompress Genital Pad Women

The Juzo SoftCompress Pad for the female genital area is worn underneath the compression garment. It provides a targeted increase in pressure, cushioning and an ideal pressure distribution when treating genital oedemas.

Description	Art. no.	Standard size	Custom-made
		1-3	
 Genital pad women size 1 – 20 x 20 cm size 2 – 30 x 28 cm size 3 – 40 x 32 cm	6839	●	●




Genital pad men

6840

Juzo SoftCompress Genital Pad Men

The reversible Juzo SoftCompress Pad for the male genital area is worn underneath the compression garment. The broad channel structure is suitable for fibrotic tissue while the narrow structure additionally promotes the lymph drainage. It provides a targeted increase in pressure and an ideal pressure distribution when treating genital oedemas.

Description	Art. no.	Standard size	Custom-made
 Genital pad men size 1 – 28 x 20 cm size 2 – 30 x 23 cm size 3 – 32 x 26 cm	6840	1-3	
		●	●




Knee joint pad

6845

Juzo SoftCompress Knee Joint Pad

The Juzo SoftCompress Knee Joint Pads cushion and soften the area in the bend of the knee joint.

Description	Art. no.	Uni size	Custom-made
 Knee joint pad 13 x 17 cm	6845		
		●	●




Elbow joint pad

6842

Juzo SoftCompress Elbow Joint Pad

The Juzo SoftCompress Elbow Joint Pads cushion and soften the area in the bend of the elbow joint.

Description	Art. no.	Uni size	Custom-made
 Elbow joint pad 10 x 13 cm	6842		
		●	●




Malleo pad

6843

Juzo SoftCompress Malleo Pad

Juzo SoftCompress Malleo Pads cushion the Bisgaard (perimalleolar) region of the ankle. They provide a targeted increase in pressure and an ideal pressure distribution. The pads are delivered in pairs.

Description	Art. no.	Uni size	Custom-made
 Malleo pad 6.5 x 9 cm delivered in pairs	6843	●	●




Hand pad

6844

Juzo SoftCompress Hand Pad

Juzo SoftCompress Hand Pads are especially for the back of the hand with two finger openings. They provide a targeted increase in pressure and an ideal pressure distribution.

Description	Art. no.	Uni size	Custom-made
 Hand pad 8.5 x 24 cm	6844	●	●

Care instructions



Material composition

Fabric: 100 % Cotton

Foam: 100 % Polyurethane

Storage instructions

Always store Juzo SoftCompress products at room temperature and away from direct sunlight.

Washing and care instructions



- Please wash Juzo SoftCompress products daily.
- Juzo SoftCompress products can be washed in the washing machine at 60 °C using a laundry detergent for delicate fabrics. Please do not use any fabric softener.
- Do not bleach Juzo SoftCompress products.
- Do not dry in a clothes dryer.
- You can reduce the drying time by placing Juzo SoftCompress products on a thick towel, rolling the towel up tightly, and firmly pressing out excess moisture. Do not leave Juzo SoftCompress products in the towel and do not dry on a radiator or in the sun.
- Do not iron Juzo SoftCompress products.
- Do not dry clean Juzo SoftCompress products.

International
Julius Zorn GmbH
Juliusplatz 1
86551 Aichach
Germany
info@juzo.com

www.juzo.com

Juzo UK Ltd
South Court
Sharston Rd
Sharston
Manchester M22 4SN
Great Britain
info@juzo.co.uk

www.juzo.co.uk

