

Product Information



**Fall prevention alarms
& anti-wandering products**



Tested & approved
To IEC 60601-1-2
Medical Devices
General



Prevent Falls & Wandering

Dear Health Care Provider,

HealthSaver have been very involved in the fall prevention/protection products field for many years with market leading products like HipSaver and HeadSaver.

Over the years we have been asked many times for good quality, well priced fall alarms and anti-wandering solutions, and in 2017 we teamed up with Smart Caregiver Corporation in the US to supply their range of alarm systems in Australia and New Zealand. All Smart Caregiver products must pass rigorous testing and exceed industry standards for quality and reliability. The Smart care mission is to "Supply the least expensive and best quality resident fall prevention monitors so that caregivers can afford the quantity of devices necessary to keep their residents safe and comfortable". A sentiment we certainly agree with at HealthSaver.

The Smart Caregiver range of fall prevention and anti-wandering products include:

- Bed Exit Monitoring
- Chair Exit Monitoring
- Door Exit Monitoring
- Motion Activated Alarms
- Anti-Wandering Solutions
- Portable Paging Systems

Quieter cordless fall prevention means no noise or tripping hazards in the resident's room. For example the Cordless® Fall Monitor can be mounted away from the resident's bed, removing in-room alarm noise and wires from the sensor to the monitor. This monitor works with bed pads, chair pads and floor mats with embedded transmitters that send a wireless signal to the fall monitor.

So whether a caregiver at home, a skilled nursing facility, or an assisted living facility needing a quiet, audible and/or visual fall prevention or anti-wandering system, Smart Caregiver Corporation and HealthSaver Pty Ltd have the solution.

As always, please contact the team at HealthSaver with questions or comments and of course technical advice and help.

Sincerely,



Ian Lancaster
Managing Director
HealthSaver Pty Ltd

Table of Contents



Smart Caregiver Alarm Monitors Features	4
Low Cost Bed & Chair Exit Alarm Monitor (TL-2100B)	5
Easy- To -Use Bed / Chair Exit Alarm Monitor (TL-2100E)	6
Safety Auto-Reset Alarm Monitor (TL-2100S)	7
Change Pad Indicator Alarm Monitor (TL-2100CP)	8
Recordable Voice Alarm Monitor Compatable With Bed/Chair Pad <u>OR</u> Pull-String (TL-3100V)	9
Quiet Fall Alarm Monitor with Wireless Signal to Pager (TL-2016R)	10
Low Cost, Quiet and CordLess & Wireless Fall Prevention Alarm Monitor (433-EC)	11
Quieter CordLess® Fall Alarm Monitor (TL-2100G)	12
CORDED Bed Pads, Chair Pads & Floor Mats For Bed & Chair Exit Alarm Monitors	13
CORDLESS® Bed Pads, Chair Pads & Floor Mats For Bed & Chair Exit Alarm Monitors	14
TIMED CORDED Bed & Chair Pads for Use With Change Pad Indicator Alarm Monitor	15
Pull-String Alarm Monitor "UNBREAKABLE" (TL-2000)	16
Non-Magnetic Pull-String Alarm Monitor (TL-01)	17
Wireless Components & Accessories	18
Wireless Multi- Channel Central Monitoring Unit (433-CMU-40)	19
Caregiver Pagers	20
Early-Warning Chair Belts (TL-2109/TL2109V)	21
Stand Alone Call Button Paging System (TL-5102TP)	22
Motion Activated Stand Alone Caregiver Paging System (TL5102MP)	23
Retail Packaged Systems	24
HipSaver® Hip Protectors	25
HeadSaver™ Head Protection	26









HealthSaver Pty Ltd, 14/140 Wecker Road, Mansfield, QLD 4122, Australia.
Phone: 1300 767 888 Fax: 1300 767 999 Email: info@healthsaver.com.au

www.healthsaveralarms.com.au



Smart Caregiver Alarm Monitors Features

	Basic Monitor (TL-2100B)	Economy (TL-2100E)	Safety Auto-Reset (TL-2100S)	Change Pad Indicator (TL-2100CP)	Dual Voice Recordable (TL-3100V)	Wireless Fall Monitor (TL-2016R)	Economy Cordless (433-EC)	Cordless (TL-2100G)	Pull String (TL-2000)	Non-Magnetic Pull-String (TL-01)
										
PAGE NO.	5	6	7	8	9	10	11	12	16	17
Corded pads	●	●	●	●	●	●		●		
Cordless pads							●	●		
Timed pads				●						
Voice recordable					●					
Adjustable tones		●	●	●	●	●		●		
Silent w/nurse call cord	●	●	●	●	●	●		●		
Adjustable volume	●	●	●	●	●	●	●	●		
Tamper resistant			● Caregiver Key	● Caregiver Key				● Caregiver Key		
Nurse call port	●	●	●	●	●	●		●		
Auto-reset	●	●	●	●	●	●	●	●		
Always-on	●		●	●	●	●		●	●	●
On-off switch		●					●			
Alert light	●	●	●	●	●	●	●	●		
Delay option			●	●		●				
Low battery indicator	●	●	●	●	●	●	●	●		
Batteries used	1 x 9V	3 x AA	3 x AA	3 x AA	3 x AA	3 x AA	3 x C	3 x AA	1 x 9V	1 x 9V
Optional AC adapter	● AC-02		● AC-05	● AC-05	● AC-05	● AC-05	● AC-04	● AC-05		
Protective case		●	●	●	●	●		●	●	
Wall mount option	●	●	●	●	●	●	●	●		
Cord strain relief	●	●	●	●	●	●		●		
Wireless Remote Reset							●			
Pull String option					●				●	●
Warranty	1 Year	1 Year	1 Year	1 Year	1 Year	1 Year	1 Year	1 Year	LIFETIME	1 Year

TL-2100B

Low Cost Bed & Chair Exit Alarm Monitor



Great Value! All The Functions

The Smart Caregiver TL-2100B Monitor is a Full-Function, Low Cost fall monitor designed for use with all Smart Caregiver Corded Pads, Floor Mats and Seat Belts. This is a high quality product at a great, everyday value.

Auto-Reset Function

Monitor can be placed at the bedside or wheelchair, connected to a corded sensor pad. The sensor pad is placed under the resident. When the resident gets up and pressure is removed from the pad, the monitor will sound alerting the caregiver. When pressure is re-applied to the pad, the monitor will silence and reset for use.

Monitor Can Be Set To Silent

Monitor can be set to silent only when using an existing nurse call system and additional nurse call cord.

Key Features:

- Full-Range Volume Control
- Nurse Call Capability - enables you to use with existing hard-wired nurse call systems
- Low Battery & Status Indicator Lights
- Battery not included (1 x 9V)
- Optional AC adapter available (AC-02)
- **Connects to one corded sensor only**

For use with only one of the following Corded Sensors:

- PPB-WI – Bed Sensor Pad 20" x 30" (50 x 76cm)
- PPB-RI – Bed Sensor Pad 10" x 30" (25 x 76cm)
- PPC-WI – Chair Sensor Pad 10" x 15" (25 x 38cm)
- TL-2109 – Easy Release Seat Belt
- TL-2109V – Quick Release Seat Belt
- FM-05 – Floor Mat 24" x 36" (61 x 91cm)
- FM-06 – Floor Mat 20" x 30" (50 x 76cm)(under rug)
- FM-07 – Floor Mat 24" x 48" (61 x 122cm)

TL-2100E

Easy-To-Use Bed / Chair Exit Alarm Monitor



Easy-to-Use Fall Monitor

This Smart Caregiver Fall Monitor Alarm functions with weight-sensing bed and chair pressure pads, floor mats, or early warning seat belts. Caregiver is notified by an audible and visual alert.

Monitor can be placed at the bedside or wheelchair and then **connected to a corded sensor pad, floor mat or seat belt**. When the resident gets up from a bed or a chair, the monitor will sound alerting the caregiver.

On/Off Switch

The fall monitor features an "On/Off" switch that prevents the alarm from sounding when the caregiver moves the resident from bed or chair, or when changing clothing or linens.

Key Features:

- Full-Range Volume Control
- Nurse Call Capability - enables you to use with existing hard-wired nurse call systems
- Low Battery & Status Indicator Lights
- Batteries not included (3 x AA)
- Protective Case
- **Connects to one corded sensor only**

For use with only one of the following Corded Sensors:

- PPB-WI – Bed Sensor Pad 20" x 30" (50 x 76cm)
- PPB-RI – Bed Sensor Pad 10" x 30" (25 x 76cm)
- PPC-WI – Chair Sensor Pad 10" x 15" (25 x 38cm)
- TL-2109 – Easy Release Seat Belt
- TL-2109V – Quick Release Seat Belt
- FM-05 – Floor Mat 24" x 36" (61 x 91cm)
- FM-06 – Floor Mat 20" x 30" (50 x 76cm)(under rug)
- FM-07 – Floor Mat 24" x 48" (61 x 122cm)

Safety Auto-Reset Alarm Monitor

TL-2100S



SAFE w/ "Always ON" Technology

The TL-2100S Monitor is designed for use with all Smart Caregiver corded sensor pads, floor mats and seat belts. Simply plug a sensor pad into the monitor and place the pad under the resident in a chair or bed. When the resident gets up and pressure is removed, the alarm will sound notifying caregivers that the resident is on the move. This alarm features a Safety Auto-Reset™ Function which will automatically reset and silence the alarm when pressure is re-applied to the sensor pad.

Optional TamperProof™ Setting

When set to TamperProof™, the Monitor alert can ONLY be reset using the provided Caregiver Key™. Selecting this setting ensures that only the designated caregiver can silence the alarm.

Key Features:

- Full-Range Volume Control
- Nurse Call Capability - enables you to use with existing hard wired nurse call system
- Optional Setting: TamperProof™ Caregiver Key allows control of who is able to reset the monitor
- Batteries not included (3 x AA)
- Optional AC adapter (AC-05)
- Protective Case
- **Connects to one corded sensor only**

For use with only one of the following Corded Sensors:

- PPB-WI – Bed Sensor Pad 20" x 30" (50 x 76cm)
- PPB-RI – Bed Sensor Pad 10" x 30" (25 x 76cm)
- PPC-WI – Chair Sensor Pad 10" x 15" (25 x 38cm)
- TL-2109 – Easy Release Seat Belt
- TL-2109V – Quick Release Seat Belt
- FM-05 – Floor Mat 24" x 36" (61 x 91cm)
- FM-06 – Floor Mat 20" x 30" (50 x 76cm)(under rug)
- FM-07 – Floor Mat 24" x 48" (61 x 122cm)

TL-2100CP

Change Pad Indicator Alarm Monitor

Eliminate The Need to Check Dates on Sensor Pads*



Keep Track Of Dates On Pads

The TL-2100CP automatically tracks pad life so caregivers no longer have to write dates on pads. The Monitor counts down the days the pad has been in use and the Check Pad Light will illuminate 10 days prior to pad expiration. A short beep will also alert caregivers that it's time to change pad. When used with Timed Bed & Chair pressure pads, this monitor will save time and staff resources!

Easy-To-Use

Monitor can be placed at the bedside or on a wheelchair and connected to a corded timed sensor pad. The sensor pad is placed under the resident and when the resident gets up and pressure is removed from the pad, the monitor will sound alerting the caregiver. Monitor will automatically reset when pressure is reapplied to the connected sensor pad.

Patented TamperProof™ Technology

This monitor has been designed with an optional TamperProof Setting which, when enabled, requires caregivers to use the provided Caregiver Key to reset the alarm. This prevents accidental tampering and ensures that only the designated caregiver can silence alerts.

Key Features:

- Full-Range Volume Control
- Nurse Call Capability - enables you to use with existing hard wired nurse call system
- Optional Setting: TamperProof™ Caregiver Key allows control of who is able to reset the monitor
- Batteries not included (3 x AA)
- Optional AC adapter available (AC-05)
- **Connects to one timed sensor only**
- Protective case

For use with one Timed Corded Sensor Pad only:

- TB-WI – Timed Bed Sensor Pad 20" x 30"
- TB-RI – Timed Bed Sensor Pad 10" x 30"
- TC-WI – Timed Chair Sensor Pad 10" x 15"

*** Must be used with Timed Sensor Pads**

TL-3100V

Recordable Voice Fall Alarm Monitor That Works With Bed/Chair Pad OR Pull-String



Unique Dual Purpose Fall Monitor

The TL-3100V DUAL Recordable Voice Alarm works as a **pull-string monitor OR sensor pad monitor** and can easily be transferred from bed to wheelchair and back. This eliminates the need for multiple monitoring devices. The alarm alerts when the magnetic pull-string is pulled free of the alarm or when pressure is removed from weight-sensing pad. Alarm resets when the magnetic pull cord is replaced, when pressure is re-applied to the connected sensor pad or when reset button is used.

Personalised Recordings

Eliminate startling alarm noise by recording a personalized message to be played back when the alarm is activated. With up to 15 seconds of playback time available, you can record a melody or familiar voice to be heard as the alert sounds.

Key Features:

- Dual Functions
- Recordable Voice Option
- Pad Lost, Low Battery & Alarm Status Indicator Lights
- Nurse Call Capability
- Monitoring Options: Corded Bed Pad, Chair Pad, Floor Mat and Seat Belt.
- Optional AC adapter available (AC-05)
- Batteries not included (3 x AA)
- Protective Case

For use with only one of the following Corded Sensors:

- PPB-WI – Bed Sensor Pad 20" x 30" (50 x 76cm)
- PPB-RI – Bed Sensor Pad 10" x 30" (25 x 76cm)
- PPC-WI – Chair Sensor Pad 10" x 15" (25 x 38cm)
- TL-2109 – Easy Release Seat Belt
- TL-2109V – Quick Release Seat Belt
- FM-05 – Floor Mat 24" x 36" (61 x 91cm)
- FM-06 – Floor Mat 20" x 30" (50 x 76cm)(under rug)
- FM-07 – Floor Mat 24" x 48" (61 x 122cm)

TL-2016R

Quiet Fall Alarm Monitor with Wireless Signal to Pager



Reduce Alarm Noise in Room!

The monitor can be set to silent to reduce in-room alarm noise. When a resident attempts to get up, a silent wireless signal is sent to a caregiver pager or central nurse's station.

Flexible & Versatile

The wireless fall monitor is designed to work with Corded Bed & Chair Sensor Pads, Floor Mats, Seat Belts, Pagers, Emergency Call Lights, Nurse Call Cords and Central Monitoring Units. This monitor allows caregivers the flexibility to monitor residents based on their unique needs.

Increase Caregiver Freedom!

Caregiver can attend to other duties and receive an alert with use of a pager when someone is getting up from a bed or chair. Can also be used with a floor mat to notify if they are approaching a door.

Key Features:

- Monitor sends wireless signal to pager up to 300 feet away (approx 100m)
- Call Button on front of monitor
- Nurse Call Capability - enables you to use with existing hard wired nurse call system
- Batteries not included (3 x AA)
- Optional AC adapter available (AC-05)
- Protective Case



TL-2016P

Caregiver Pager
 (Sold Separately)

For use with only one of the following Corded Sensors:

- PPB-WI – Bed Sensor Pad 20" x 30" (50 x 76cm)
- PPB-RI – Bed Sensor Pad 10" x 30" (25 x 76cm)
- PPC-WI – Chair Sensor Pad 10" x 15" (25 x 38cm)
- TL-2109 – Easy Release Seat Belt
- TL-2109V – Quick Release Seat Belt
- FM-05 – Floor Mat 24" x 36" (61 x 91cm)
- FM-06 – Floor Mat 20" x 30" (50 x 76cm)(under rug)
- FM-07 – Floor Mat 24" x 48" (61 x 122cm)

433-EC

Low Cost, Quiet, CordLess & Wireless Fall Prevention Alarm Monitor



No Cords! Place monitor where you need it!

Low Cost Solution

The Economy CordLess® Exit Alarm provides affordable and effective cordfree monitoring with a range up to 300 ft (approx. 100m). Ideal for home or small facility use.

Easy To Use

The Economy CordLess® Exit Alarm is easy-to-use and works with all CordLess® Sensor Pads & Floor Mats, Nurse Call Buttons, and Motion Sensors. Monitor is light-weight and can be paired with up to six (6) wireless components to serve as a portable caregiver alert.

Remove Noise From the Room

Monitor can be mounted away from the bed, removing in-room alarm noise. Adjustable volume and gentle chime alert won't startle residents and offers QUIETER fall prevention.

Key Features:

- Affordable & effective wireless monitoring
- Easy set up and operation
- **Monitor up to six (6) components with one alarm**
- Full Range Volume Control
- Batteries not included (3 x C)
- Optional AC adapter available (AC-04)

For use with Cordless Sensors (up to 6):

- 433-NC Nurse Call Button
- 433-RB Remote Reset Button
- 433-MS Motion Sensor
- 433-EXT Window/Door Exit Alarm
- GBT-WI CordLess® Bed Pad 20" x 30" (50 x 76cm)
- GBT-RI CordLess® Bed Pad 10" x 30" (25 x 76cm)
- GCT-WI CordLess® Chair Pad 10" x 15" (25 x 38cm)
- FMT-05C CordLess® Floor Mat 24" x 36" (61 x 91cm)
- FMT-06C CordLess® Floor Mat 20" x 30" (50 x 76cm)
- FMT-07C CordLess® Floor Mat 24" x 48" (61 x 122cm)

TL-2100G

Quieter CordLess® Fall Alarm Monitor



Eliminates In-Room Noise

The CordLess® Fall Monitor is cord free and can be mounted where you need it - away from residents - eliminating in-room alarm noise. Monitor can be easily wall-mounted at the bedside or in a hallway up to 30 ft (10m) away using the provided wall-mounting hardware.

Patented Safety Features

This monitor features Patented, cutting-edge technology that automatically lets a caregiver know if the signal between pad and monitor is lost or if the monitor is out of range. In addition, this monitor features an optional Patented TamperProof™ setting, which when enabled requires caregivers to use the provided Caregiver Key™ to reset the alarm. This ensures that the monitor cannot accidentally be reset.

No Cords!

Cordless alarm system has NO cords between pads and monitors, which reduces tripping hazards and eliminates broken or tangled pad cords.

For use with up to 2 Cordless Sensor Pads:

- | | |
|--|--|
| • GBT-WI – Cordless Bed Pad 20" x 30"
(50 x 76cm) | • FMT-05C – Cordless Floor Mat 24" x 36"
(61 x 91cm) |
| • GBT-RI – Cordless Bed Pad 10" x 30"
(25 x 76cm) | • FMT-07C – Cordless Floor Mat 24" x 48"
(61 x 122cm) (Under rug) |
| • GCT-WI – Cordless Chair Pad 20" x 30"
(50 x 76cm) | • FMT-06C – Cordless Floor Mat 20" x 30"
(50 x 76cm) |

Key Features:

- **QUIETER** – Fully adjustable volume. Place monitor away from resident (i.e. above bed or in hallway)
- No cords between pads and monitors - reduce tripping hazards and eliminates broken or tangled pad cords
- Patented Pad Signal Lost notification - alerts caregiver if sensor pad and monitor are no longer communicating
- Optional Setting: TamperProof™ Caregiver Key allows control of who is able to reset the monitor
- Raised Visual Status light enabling caregiver to see alarm from a distance
- Monitor up to two components at one time
- Batteries not included (3 x AA)
- Optional AC adapter available (AC-05)
- Protective Case

CORDED Bed Pads, Chair Pads & Floor Mats For Bed & Chair Exit Alarm Monitors



Simple early warning system to help prevent falls.

Unique Safe T-Release cord separates instead of breaking or pulling away from the pad.

Bed and chair sensor pads alert caregivers whenever pressure is removed from the sensor pad (i.e. when a resident gets up from a bed or chair and pressure is removed from the pad).

Floor Mats alert caregivers when pressure is applied to the sensor mat (i.e. when a resident steps onto a floor mat and pressure is placed on the mat).

Floor Mats are ideal for bedside and doorway mobility monitoring. Bed and chair pads are ideal early warning alerts of resident mobility as they exit beds and chairs.

MFR#	SMART PRODUCT DESCRIPTION
PPB-WI	Corded Bed Pad 20" x 30" (50 x 76cm) (1 year pad life)
PPB-RI	Corded Bed Pad 10" x 30" (25 x 76cm) (1 year pad life)
PPC-WI	Corded Chair Pad 10" x 15" (25 x 38cm) (1 year pad life) (No Safe T-Release Cord)
FM-05	Corded Floor Mat 24" x 36" Grey (61 x 91cm) (1 year warranty)
FM-06	Corded Floor Mat 20" x 30" Grey for under rug use (50 x 76cm) (1 year warranty)
FM-07	Corded Floor Mat 24" x 48" Grey (61 x 122cm) (1 year warranty)

Works with: Corded Smart Fall Alarm Monitors



TL-2100B
Basic Monitor



TL-2100E
Economy



TL-2100S
Safety Auto-Reset



TL-3100V
Dual Voice Recordable



TL-2016R
Wireless Fall Monitor

CORDLESS® Bed Pads, Chair Pads & Floor Mats For Bed & Chair Exit Alarm Monitors



No Cords Between Pads & Monitors

When using a CordLess® Pad with a CordLess® Monitor you are able to remove the alarm noise from the room, creating a quieter environment for residents.

The Pad is placed under resident, either on a chair or bed, and the floor mat can be placed either by the bed or by the door. When a resident gets up from a pad, or steps on a floor mat, a wireless signal is sent to the monitor setting off an alarm to alert caregivers.

MFR#	SMART PRODUCT DESCRIPTION
GBT-WI	CordLess® Bed Pad 20" x 30" (50 x 76cm) (1 year pad life)
GBT-RI	CordLess® Bed Pad 10" x 30" (25 x 76cm) (1 year pad life)
GCT-WI	CordLess® Chair Pad 10" x 15" (25 x 38cm) (1 year pad life)
FMT-05C	CordLess® Floor Mat 24" x 36" Grey (61 x 91cm) (1 year warranty)
FMT-06C	CordLess® Floor Mat 20" x 30" Grey for under rug use (50 x 76cm) (1 year warranty)
FMT-07C	CordLess® Floor Mat 24" x 48" Grey (61 x 122cm) (1 year warranty)

Works with: Cordless Smart Fall Alarm Monitors



433-CMU
 Central
 Monitoring Unit



433-EC
 Cordless®
 Alarm



TL-2100G
 Cordless®
 Monitor

Timed Corded Bed & Chair Pads for Use With Change Pad Indicator Alarm Monitor



When used with the Change Pad Indicator Monitor, the Monitor automatically tracks the life of the sensor pad and alerts caregiver 10 days prior to Pad expiration date.

Timed Chair Pads

- Unique coiled cord reduces tripping and entanglement hazards

Timed Bed Pads

- Unique Safe-T-Release cord separates instead of breaking or pulling away from the pad

MFR#	SMART PRODUCT DESCRIPTION
TB-WI	Timed Bed Pad 20"x30" (50 x 76cm) (1 year warranty)
TB-RI	Timed Bed Pad 10"x30" (25 x 76cm) (1 year warranty)
TC-WI	Timed Chair Pad 10"x15" (25 x 38cm) (1 year warranty)

Works only with: Change Pad Indicator Fall Alarm Monitor



TL-2100CP
Change Pad
Indicator Fall
Monitor

TL-2000

Lifetime Pull-String Alarm Monitor. If it breaks, we replace it!



Lifetime Warranty

No more replacing expensive pull-string monitors! Our policy is simple: if it breaks, we replace it - FREE!


The virtually indestructible protective cover will stop the monitor from breaking when dropped. Many facilities still like the simplicity of a pullstring but hate having to replace them all the time.

Locking Mounting Strap

The integrated, locking mounting strap easily and securely mounts to chair backs. Clip attaches to clients clothing. Alarm activates upon separation.

Key Features:

- Unbreakable
- Affordable
- Magnet alarm resets when magnet is replaced
- Thick, impact resistant foam cover
- Battery included (1 x 9V)
- Unbreakable Protective Casing
- Stand-alone unit

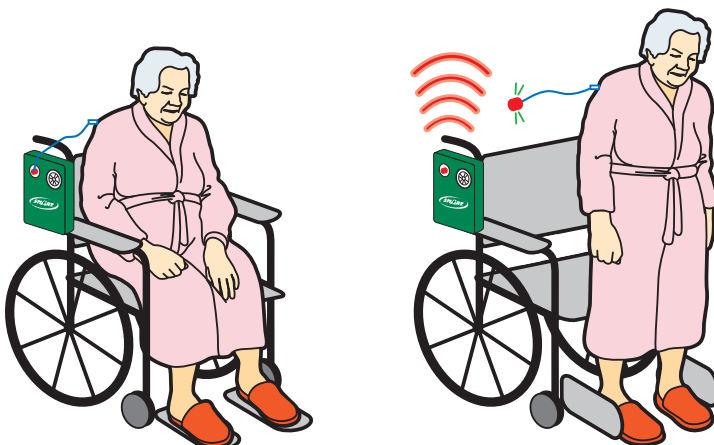

Non-Magnetic Pull-String Alarm Monitor

Pull-string without the magnetic interference that could affect pacemakers

TL-01



How It Works


This Pull-String Fall Monitor features clip-positioned cord that attaches via a garment clip to the resident's clothing. When the resident attempts to get out of their chair, the tension on the Pull-String Cord causes the clips to pull away from their position, causing the fall alarm to sound, alerting the caregiver of the resident's activity. Adjust the Cord Stop to desired Pull-String length for the resident's comfort and also to prevent false alarms.

Ideal for residents with electronic implants


This is perfect for any residents with a pacemaker. It will still provide the safety of a traditional pull-string without the magnetic interference that could affect pacemakers.

Key Features:

- No Magnets. Won't interfere with magnet-sensitive equipment
- A breakaway buckle is included for easy mounting
- Battery included (1 x 9V)
- Stand-alone unit
- 1 Year Warranty

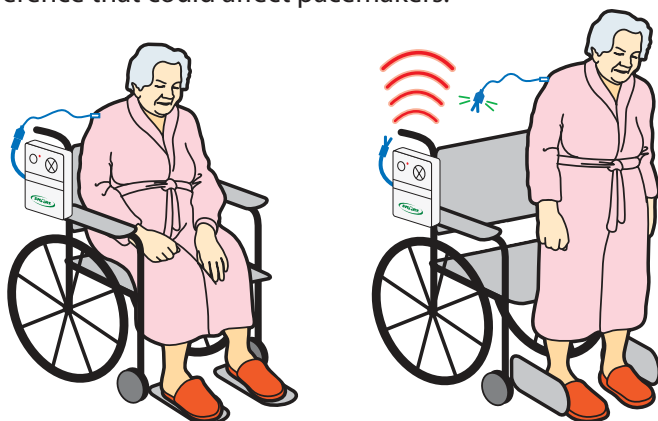

Wireless Components & Accessories

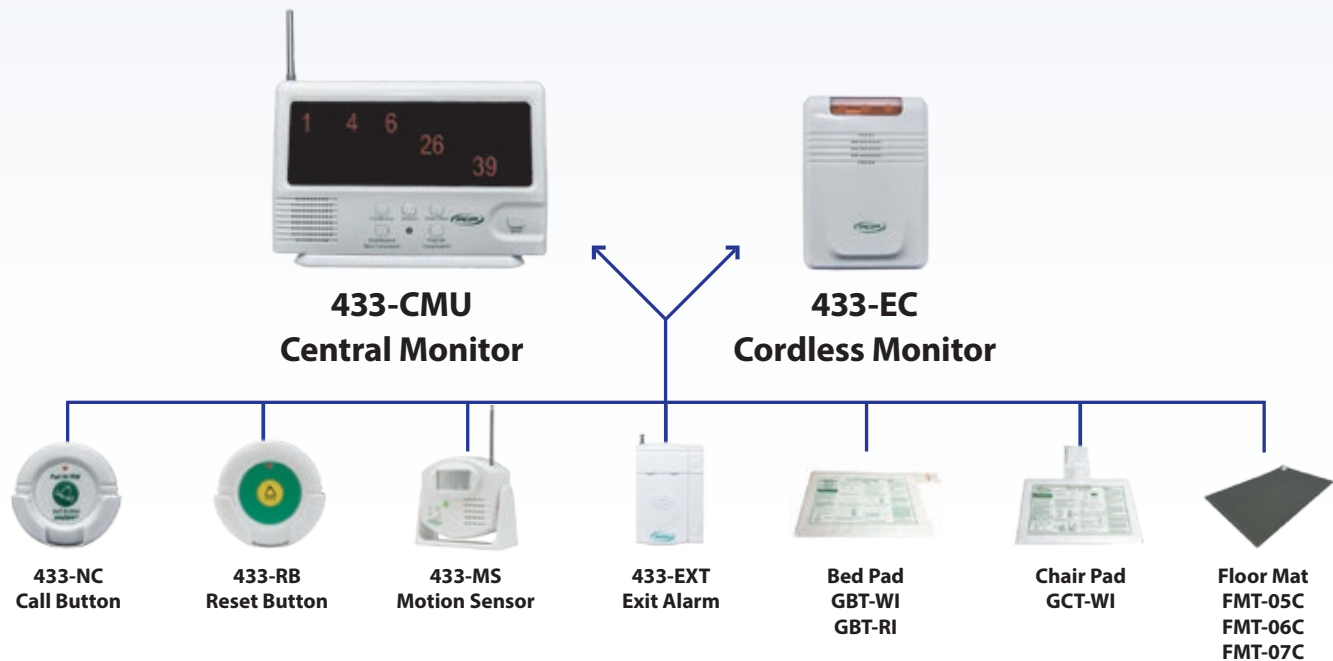

433-NC Nurse Call Button

The 433-NC Wireless Nurse Call Button allows individuals to signal caregivers for assistance with push-button activation. Call button comes standard with necklace lanyard and wall-mounting cradle.

433-RB Remote Reset Button

The 433-RB Wireless Remote Reset Button allows staff to reset the alarm monitor remotely while attending directly to resident needs.

433-MS Motion Sensor


The 433-MS Wireless Motion Sensor has flexible mounting options to capture resident movement from any angle. Optional AC Power Adapter available.

433-EXT Window/Door Exit Alarm

The 433-EXT Wireless Window/Door Exit Alarm easily mounts to window or door and sends a signal to monitor when window/door is opened and the connection is broken.

433-PGD/PRB Caregiver Pager

The 433-PGD/PRB Caregiver Pager is pocket sized and designed for use with the 433 Central Monitoring Unit. The CMU sends a signal to the Caregiver Pager and the large LCD screen clearly indicates which resident/device triggered the alarm. 433-PRB has a reset button that resets the 433-CMU-40/60 when alarming.

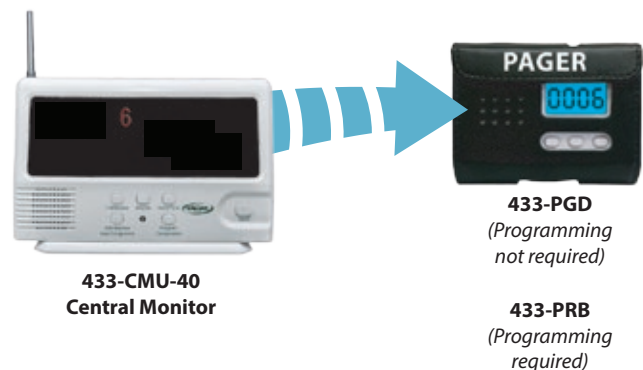


433-CMU-40



Wireless Multi-Channel Central Monitoring Unit

Monitor multiple residents from one central location

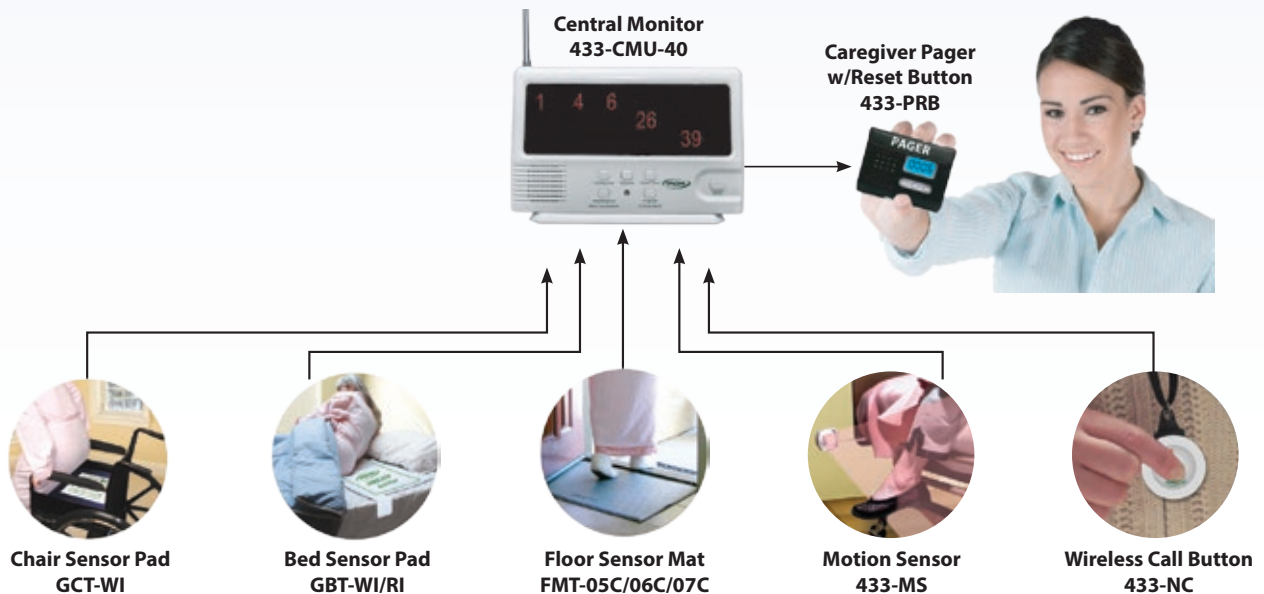

Eliminate In-Room Alarm Noise

The Smart Caregiver 433-CMU Central Monitoring Unit allows you to monitor multiple wireless components from one central location, such as a Nurses' Station. This eliminates in-room alarm noise & significantly reduces the tripping hazards associated with corded units.

Know Who Needs Help

The Central Monitor alerts audibly and the display visually indicates which resident needs assistance. At a glance, this provides caregivers with the information they need to help to protect their residents from falls.

Local Pager Notification

The CMU communicates with any number of 433-Pager Units, and upon activation, the pager audibly alerts and the LCD screen displays which component activated the alert. This enables caregivers to attend to their duties yet be alerted in real time that a resident needs assistance.

433-PRB has a reset button that resets the 433-CMU-40 when alarming.

For use with Caregiver Pagers:

- 433-PGD (Programming not required)
- 433-PRB (Programming required)

Key Features:

- Alarms audibly & visually at a central location and portable caregiver pager(s)
- Compatible with all Smart wireless components
- 40 & 60 channel capacity models available
- Range up to 300 feet (approx. 100m)
- Systems may be pre-programmed prior to delivery or easily re-programmed using the provided instructions
- Optional AC adapter available (AC-04)
- Batteries not included (4 x C)

Caregiver Pagers



Central Monitor
433-CMU-40



Caregiver Pager
433-PGD
433-PRB



Caregiver Pager
TL-2016P



Wireless Monitor
Unit
TL-2016R

433-PGD

Install 2 "AA" batteries. Turn pager on - slide power switch to upper position for chime alert or lower position for vibrate alert.



Volume
Voice/Tone
Review

The 433-PGD will work automatically with any 433-CMU Central Monitor Unit within range (up to 70m). The LCD display on the pager will show, by number, which device triggered the alert.

No programming required.



Supplied with case

433-PRB

Install 2 "AA" batteries. Turn pager on - slide power switch to upper position for chime alert or lower position for vibrate alert.



Volume
Reset
Review

The 433-PRB will work with any 433-CMU Central Monitor Unit, when programmed into that unit, within range (up to 70m). The LCD display on the pager will show, by number, which device triggered the alert.

Must be programmed into the 433-CMU.
(Refer 433-CMU instructions)



Supplied with case

TL-2016P

Install 2 "AA" batteries. Turn pager on - slide power switch to upper position for chime alert or lower position for vibrate alert.



Volume
Voice/Tone
Review

The TL-2016P will work automatically with any TL-2016R Monitor Unit within range (up to 70m). The LCD display on the pager will show, by number, which device triggered the alert.

No programming required.



Supplied with case

TL-2109/TL2109V

Early-Warning Chair Belts

Know before resident gets up!



For Use With Wheelchair Exit Alarms

These early warning seat belts are designed for use with wheelchairs and all of our Corded Fall Monitors.

Early Warning

The Fall Prevention Wheelchair Seat Belts reduce falls by triggering an alert when the Easy Release Buckle is unbuckled or the Quick Release Hook & Loop strap is opened. This notifies a caregiver before a resident attempts to leave their chair.

Non-Restraining

Both the Easy Release and Quick Release seat belts preserve resident's freedom of movement by not acting as a restraint. The Quick Release Seat Belt is ideal for residents with low dexterity as the Hook & Loop strap can be easily unfastened with one hand.

IMPORTANT: Seat Belts are not designed to restrain or hold individuals in position in their chairs.

Key Features:

- Easy Release Seat Belt features a large, red button that is easy to locate and release with one finger
- Quick Release Seat Belt features an easily grasped Hook & Loop strap which unfastens quickly even for those with limited dexterity
- Easy Installation



**Easy Release
TL-2109**
Adjustable
20" to 54"
(51 to 137cm)



**Quick Release
TL-2109V**
Adjustable
20" to 44"
(51 to 111cm)

For use with one of the following Alarm Monitors:



TL-2100B
Basic Monitor



TL-2100E
Economy



TL-2100S
Safety Auto-Reset



TL-3100V
Dual Voice Recordable



TL-2016R
Wireless Fall Monitor

TL-5102TP



Stand Alone Call Button Paging System



Easy-To-Use Two Call Button Paging System!

Increased Caregiver Freedom

The TL-5102TP is a Stand Alone Personal Paging System and includes two Call Buttons and one Caregiver Pager. When pressed, the Nurse Call Buttons send a wireless signal to the pager, alerting the caregiver that assistance is required. This simple system eliminates in-room alarm noise and with a range of up to 300 feet (approx. 100m), gives caregivers the freedom they need.

Easy Push-Button Operation

The call buttons can be worn as a necklace pendant using the provided lanyard strap or wall-mounted using the provided cradle for easy access.

Versatile Pocket-Sized Pager

The Caregiver Pager has two audible sounds to choose from and is small enough to be carried in a pocket or clipped to a belt.

Key Features:

- Insert batteries and it's ready to use
- Audible alert pager with optional tones
- Range up to 300 feet (approx. 100m)
- Lanyard Strap & Wall-Mounting Cradle included
- Batteries included (2 x AA, 2 x 12V)

"Best thing ever!"

It enables me to do other stuff and my 80 years old mother can page me when she needs my assistance."

- Colleen B.
Palm Coast, FL.

TL-5102MP

Motion Activated Stand Alone Caregiver Paging System



Easy-To-Use Caregiver Alert!

No Noise in Room!

The Smart Caregiver TL-5102MP is a stand-alone motion paging system that allows you to eliminate in-room alarm noise. The Motion Sensor is placed in a doorway or by the bedside. When motion is detected, a wireless signal is sent to the caregiver pager up to 300 feet (approx. 100m) away.

Caregiver Alerted by Pager

The Caregiver Pager is pocket-sized and light-weight; giving caregivers the freedom they need.

Place it where you need it

The sensor can be placed next to the bed, by a doorway or anywhere to detect movement.

Key Features:

- Quiet - No alarming monitors
- Easy to set up - insert batteries and it's ready to use
- Pocket sized pager - Easy to carry
- Motion Sensor can be placed by bed or doorway to monitor residents activities
- Batteries included (2 x AA, 3 x AA)
- Optional AC adapter available (AC-04)

Retail Packaged Systems



Two Call Buttons w/Pager
 RP-TL-5102TP



Motion Sensor to Pager
 RP-TL-5102MP



Stand-Alone Motion Sensor Alarm
 RP-TL-2700



Complete CordLess® Bed Exit Monitoring System
 RP-433BR1-SYS



Complete CordLess® Chair Exit Monitoring System
 RP-433BR1-SYS



Complete Bed Exit Monitoring System
 RP-BBR1-SYS



Complete Chair Exit Monitoring System
 RP-BC1-SYS

Also available...



HipSaver® Hip Protectors

Protects vulnerable hips during a fall.

- Soft, sewn-in, tapered pads
- Proven in biomechanical testing and clinical studies
- Machine wash & dry even in hot institution laundries



Nursing Home™

Wide leg holes, roomy fit, easy to pull up and down.



QuickChange™

QuickChange has a drop down panel with only 2 snaps for easy changing of incontinence products.



SlimFit™

Sleek design - discreet under clothing, ideal for fashion conscious residents and outpatients




Open-Bottom™

Open crotch design for easy self-toileting. Use alone or residents can wear their own underwear over it.



Tailbone Pad

Wide leg holes, roomy fit, easy to pull up and down.

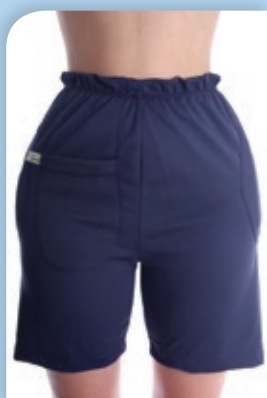


Shorts

Proven HipSaver protection in soft comfortable trackpants and shorts.

- Six hip sizes and 2 lengths
- No elastic or zippers at ankle
- Trackpants & shorts - Navy colour
- Interims - Grey colour.



Interim

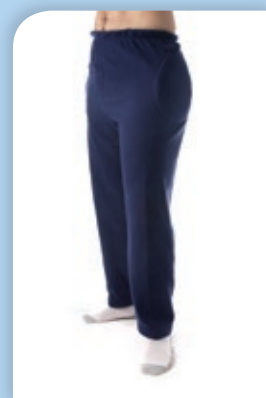

Track Pants

Available with:

- Hip Pads
- Hip & Tailbone pad
- Hip & Knee pads
- Hip, Tailbone & Knee pads

Also available...



HeadSaver® - The Right Choice for Head Protection

Comfort and ease of use are just two of the benefits enjoyed by all **HeadSaver** users, from Nursing Home residents and staff to active, independent seniors.

Although we emphasize the importance of our head protection products for the elderly, **HeadSaver** is the right head protector choice for all with concerns about falling or with medical conditions such as Epilepsy, Parkinson's Disease, or Cerebral Palsy, which make them unsteady on their feet.

HeadSavers are available for both male and females, designed to fit the lifestyles and specific needs of those who use them. They are available in two sizes which are fully adjustable to fit most people. When placed on the head, HeadSaver offers an extremely comfortable fit that holds the HeadSaver in place so that during a fall the HeadSaver should not come off the head. If additional securing to the wearer's head is required a detachable chin strap with clip is provided with every HeadSaver.



Available in 2 sizes;
Small/Medium & Medium/Large.

Brain injury is irreversible!

Don't take the risk.

Reduce the risk of devastating head and brain injury with HeadSaver - the stylish, lightweight, comfortable, all-soft head protector.



So why shouldn't a head protector look like any ordinary headwear?

We believe that dignity (and fashion!) are everyone's right. That's why we have designed and developed various head coverings to "dress up" the HeadSaver soft head protector. Whether it's a sun hat, winter beanie or a fashionable scarf, they are all available as optional items, are easily fitted and can be changed in an instant.



Scarf (Blue or Pink)



Sun Hat



Beanie





Important Notes

Please read the instructions and the Legal Product Disclaimer supplied with each product before using each Smart Caregiver device.

Smart Caregiver products are designed to be installed by the end user. As such it is, entirely the responsibility of the buyer to ensure the products are properly installed and tested.

Further, Smart Caregiver products are not intended to replace good caregiver practice including, but not limited to, direct patient supervision, staff training and testing of the system before each use.

WARNING: Smart Caregiver and/or Healthsaver Pty Ltd does not claim that the devices will stop elopement and/or stop falls. The devices are designed to augment caregiver's comprehensive resident mobility management programme.

Notes



*Fall prevention alarms
& anti-wandering products*



HealthSaver Pty Ltd, 14/140 Wecker Road, Mansfield, QLD 4122, Australia.
Phone: 1300 767 888 Fax: 1300 767 999 Email: info@healthsaver.com.au

www.healthsaver.com.au

