




hygiena

Helping Hospitals Define Clean

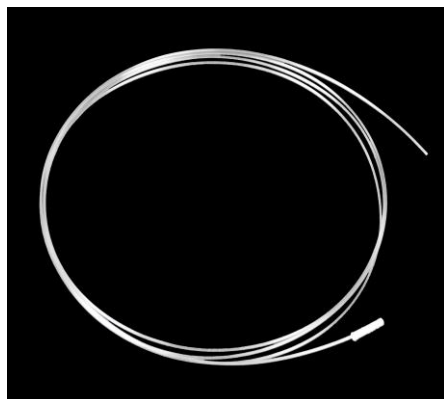


EndoSwab | Endoscope Swab

Hygiena EndoSwab enables cleaning verification of scopes or other reusable instruments with hard to reach internal channels. The device features an absorbent foam tip attached to a long flexible wand. After sampling, the swab tip is removed and used with Hygiena ATP tests to detect any residual contamination left behind after cleaning, ensuring scopes and other reusable instruments are ready for disinfection and cleaning procedures were effective.



Find support documents, instructional videos, and more at www.hygiena.com

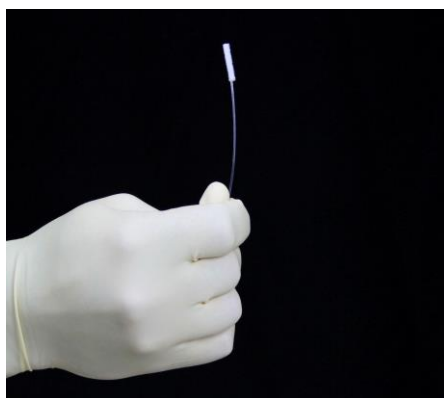


Benefits:

- Results in 15 seconds
- 3 mm absorbent sponge tip for maximum sample recovery
- 2.4 m wand length ensures recovery of a representative sample for most endoscope lengths
- Clinically clean packaging
- Individually bagged for single use

Compatibility:

- Use with Hygiena UltraSnap ATP tests or SuperSnap high sensitivity ATP tests, depending on desired sensitivity

EndoSwab Instructions:

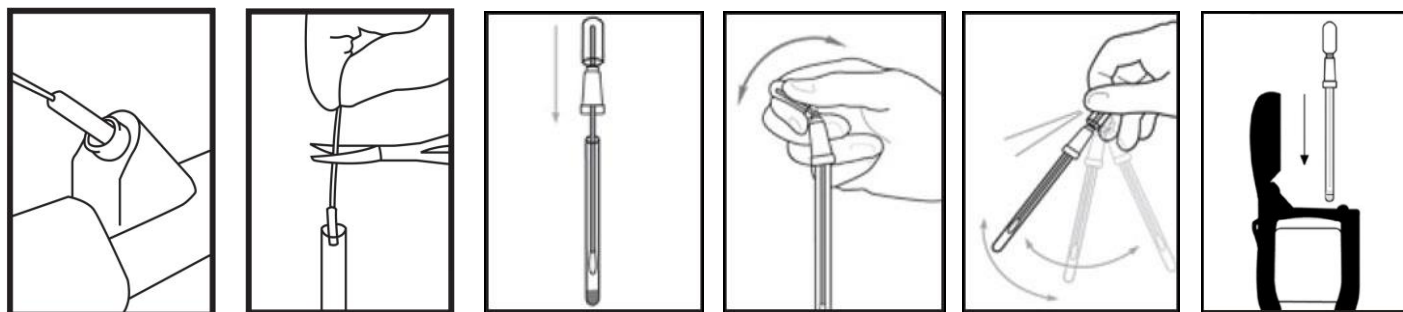

Catalog No.	Description	Quantity
ENDOSWAB	Endoswab – Swab for Endoscope Cleaning Verification	100
US2020	UltraSnap ATP Swabs	100
SS3H	SystemSURE Plus ATP Cleaning Verification System	1



ANAERON
PTY LTD

T: 02 8209 3815 E: info@anaeron.com.au
www.anaeron.com.au

Americas Tel: +1.805.388.8007 | Int

www.hygiena.com | Email: info@hygiena.com

REV C 052015