

# eue

Ace of Face

*Hydrodermabrasion, Diamond Peeling,  
Microcurrent & EMS Device*



 TADLI The Australian Dermal  
& Laser Institute

# Features & Benefits



## Variety of Microcurrent Applicators

Ball shape , T-shape and comb-shape microcurrent metal tips, for convenience in various treatment areas.

Conductive micro-gloves, for a more pleasant and attentive microcurrent treatments.

## Multiple Hydro/Diamond Tips



3 different diamond tips, in different densities and sizes, for exfoliation of different skin types and facial areas.

## Control Dials

Adjust the hydrodermabrasion jet of X-Stream toner and the vacuum suction, per patients' skin type, condition and individual preferences.



## Hydrodermabrasion Toner

The exclusive formula of X-Stream, the Eve hydro-dermabrasion toner, is highly effective in **gently exfoliating and effectively hydrating** the skin.



## Smart Toner / Sediment Mechanisms

Smart mechanism enables **easy plugging/un-plugging** the bottle of X-Stream toner.

Smart sensor presents on-screen warning, to **prevent overflow of the sediment bottle**.

## Structured, yet Adjustable, Treatment Programs

A selection of **pre-set, automated programs**, while allowing **adjustment of program duration and parameters** to the treatment indications, skin conditions, and patients' preferences.

# Hydrodermabrasion


Hydro dermabrasion by Eve thoroughly cleans, gently exfoliates and well-hydrates the skin surface. It is the ideal first step in every aesthetic treatment.

The result of this treatment is a healthier, well-hydrated, fresh looking and glowing skin, with a much smoother texture. You will also notice significant reduction of blackheads and enlarged pores.



## Cleansing

The X-Stream jet removes dirt, excess sebum and dead skin cells from the skin surface.

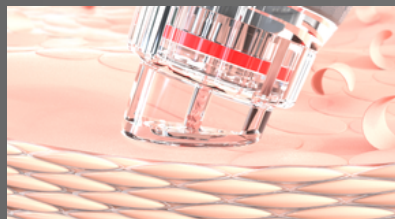

## Exfoliation

PHA, LHA, and proteolytic enzymes enable very gentle exfoliation of dead skin cells and waste.



## Extraction

Adjustable vacuum suction painlessly extracts impurities from inside the pores.



## Hydration

The ingredients of the X-Stream toner effectively hydrate the stratum corneum.



## Eve X-Stream Hydradermabrasion Solution

- **PHA (Poly Hydroxy Acid) – Lactobionic Acid:** Effective moisturiser and mild exfoliator, even for sensitive skin. Proven to effectively improve texture and complexion.
- **LHA (Lipo Hydroxy Acid) – Capryloyl Salicylic Acid:** Enables gentle, non-irritating exfoliation, while reducing redness and rejuvenating the skin surface. Especially effective in decreasing P. acne bacteria in the sebaceous glands.
- **Proteolytic Enzymes – Bromelain:** Helps the stratum corneum's own proteolytic enzymes to break down protein Corneodesmosomes. This process is essential for proper desquamation and cell turnover in the stratum corneum.
- **Hyaluronic Acid & Niacinamide:** Improves skin texture, revitalises dull complexion, lightens pigmentation, restores the natural water/oil balance, soothes irritation, anti-oxidant effect.

# *Diamond Peeling Microdermabrasion*

Effortlessly eliminate dead skin cells, impurities and excess sebum, revealing a smoother, brighter complexion. With adjustable vacuum suction and three interchangeable diamond tips, this function can be used as a stand-alone treatment or as a pre-treatment for other services.

There are a selection of hydrodermabrasion tips, for different facial areas and for the scalp.



## *Electrical Muscle Stimulation (EMS)*

A "training session" for the facial muscles, as EVE delivers low-intensity, low-frequency electric currents into the skin, causing repeated muscle contraction. These microcurrents mimic the body's natural bioelectric currents, regulating and reigniting various intracellular electro-dependent processes.

Resulting in reactivated muscle function, restored muscle tone, and a visible lifting and tightening effect that enhances facial contours and overall appearance.









# Microcurrent

The Eve microcurrent function delivers low intensity and low frequency electric impulses into the skin. Though hardly perceptible, these microcurrents mimic the body's own bioelectric currents, and hence regulate/re-ignite electro-dependent intracellular processes.


Eve offers 4 microcurrent functions, each differing in intensity and frequency of electric currents. Each microcurrent function regulates the bioelectric activity of different target cells, therefore having a different physiological impact on the skin:

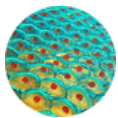

**MTB:** Re-activating Skin Cell Metabolism

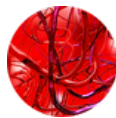

**LMP:** Stimulating Micro-Circulation & Lymphatic Drainage

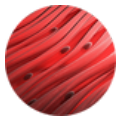

**MYO:** Myo-Stimulation & Myo-Relaxation



**ALV:** Alleviating Skin Pain & Discomfort

## Variety of Microcurrent Applicators

- **Ball Shape:** Suitable for any microcurrent function and for all treatment areas.
- **Comb Shape:** Intended for the scalp, to stimulate micro-circulation, thus strengthen the hair roots and reduce seborrhea and dandruff.
- **T-Shape:** Useful for covering broader skin areas (e.g. the cheeks or the neck) conveniently and more efficiently.
- **Conductive Gloves:** For a pleasant, attentive and highly-effective microcurrent treatment.



# About TADLI

*Raising the standard of skin education & technology and bridging the gap between beauty, dermal & health care practitioners, ensuring an integrated approach.*



Established in 2007, The Australian Dermal & Laser Institute (TADLI) offer diverse and comprehensive aesthetic device training for various manufacturers and direct to clinics. We believe that "best practice" is not negotiable and safe, superior results are the minimum standard. This is key to having credible staff members, ensuring greater client satisfaction, and guaranteeing a successful business.

TADLI are proud to offer an array of aesthetic devices that are both cost effective and results driven. Our strong focus on education and maximising treatment results, makes us a trusted partner to ensure clinics receive a comprehensive experience.



#### **Award-Winning Training:**

Every device includes award-winning training, providing users the expertise to maximise potential.



#### **Continuous Clinical Support:**

We stand by our partners with ongoing clinical support, providing insights, solutions and advice whenever required.



#### **Trust and Integrity:**

Our core values emphasise building enduring relationships based on trust, making us a dependable partner in your clinic's aesthetic journey.



#### **Promotion of Safe Practices:**

With our extensive training and continuous support, clinics can ensure they're always delivering treatments that are safe and compliant.



#### **Authentic Female Leadership:**

Founded and run by women who've thrived in their own aesthetic clinics, we prioritise genuine care and industry expertise over mere profit-making.

# Specifications

Dimensions (Net)	Device: 20 x 40 x 30 cm Cart: 36 x 40 x 105 cm
Dimensions (Packaging)	Device: 35 x 55 x 50 cm Cart: 48 x 18 x 76 cm
Weight (Packaging)	Device: 13kg Cart: 8.5kg
Device Power Consumption	100-240V, 50 / 60 Hz Max 100VA (24V / 5A)
Microcurrent Programs	Electric Current: 10-1000 $\mu$ A Frequency: 6-500 Hz Pulse Amplitude: Max 1000 $\mu$ A
EMS Programs	Electric Current: 1-3mA Frequency: 2-2000 Hz Pulse Amplitude: Max 3 mA
Hydrodermabrasion Suction Power	650 mmHg



Distributed by The Australian Dermal & Laser Institute (TADLI)  
info@dermalandlaser.com.au | (03) 8001 5004 | www.dermalandlaser.com.au